

Bring joy and happiness
to your loved ones



- Lowerlimbs prosthetics components
- Suspension system
- Cosmetics, Tools and Materials
- Orthotics and Rehabilitation products
- Upperlimbs prosthetics components

Contents

<i>Index</i>	P.3-8
<i>Regal Pediatric</i>	P.9-16
<i>Regal Feet and Adaptors</i>	P.17-32
<i>Regal Kneejoints and Hip Joints</i>	P.33-50
<i>Regal Suspension Systems</i>	P.51-68
<i>Regal Cosmetic, Tools and Materials</i>	P.69-92
<i>Regal Orthotics and Rehabilitation</i>	P.93-102
<i>Regal Upperlimbs Modular Prosthesis</i>	P.103-105
<i>Regal FAQ & Warranty</i>	P.106-107

Our Quality Management System and Product Certification:



ISO13485:2003
Design and manufacture of
Silicone Cosmetic Prostheses
and Supply of Prosthetics and
Orthotics products



ISO9001:2008
Quality management system
certification



CE certification

Index

Product Code Pediatric	Page no.
0-00-4S30C	12,36
0-00-4S40C	14,35
0-00-4S40C=LA	14,35
0-00-4S40C-01	14,35
0-00-4S40C-02	14,35
0-00-4S40C=LA-01	14,29,35
0-00-4S41C	13,34
0-00-4S41C=LA	13,34
0-00-4S41C-01	13,34
0-00-4S41C=LA-01	13,29,34
0-00-H1C	15,36
0-00-H1C-01	15,36
0-01-FAC	15,70
0-01-FBC	16,71
0-01-FTT-L-12	10,19
0-01-FTT-R-12	10,19
0-01-FTT-L-13	10,19
0-01-FTT-R-13	10,19
0-01-FTT-L-14	10,19
0-01-FTT-R-14	10,19
0-01-FTT-L-15	10,19
0-01-FTT-R-15	10,19
0-01-FTT-L-16	10,19
0-01-FTT-R-16	10,19
0-01-FTT-L-17	10,19
0-01-FTT-R-17	10,19
0-01-RLOCC-AL	16,99
0-01-RLOCC-SS	16,99
0-01-OFFSC-SS	16,99
0-02-2N1SC	12,23
0-02-2N1TC	12,23
0-02-A10AC	11,22
0-02-A1AC	10,22
0-02-A1SC	10,22
0-02-A22C	11,23
0-02-A5AC	11,22
0-02-A6AC	11,22
0-02-BC150	12,25
0-02-BC300	12,25

Product Code Feet & Adaptors	Page no.
1-01-DYNA-L-22	21
1-01-DYNA-R-22	21
1-01-DYNA-L-23	21
1-01-DYNA-R-23	21
1-01-DYNA-L-24	21
1-01-DYNA-R-24	21
1-01-DYNA-L-25	21
1-01-DYNA-R-25	21
1-01-DYNA-L-26	21
1-01-DYNA-R-26	21
1-01-DYNA-L-27	21
1-01-DYNA-R-27	21
1-01-DYNA-L-28	21
1-01-DYNA-R-28	21
1-01-DYNA-L-29	21
1-01-DYNA-R-29	21
1-01-DYNA-01	21
1-01-FT-L-21	20
1-01-FT-R-21	20
1-01-FT-L-22	20
1-01-FT-R-22	20
1-01-FT-L-23	20
1-01-FT-R-23	20
1-01-FT-L-24	20
1-01-FT-R-24	20
1-01-FT-L-25	20
1-01-FT-R-25	20
1-01-FT-L-26	20
1-01-FT-R-26	20
1-01-FT-L-27	20
1-01-FT-R-27	20
1-01-FT-L-28	20
1-01-FT-R-28	20
1-01-FT-L-29	20
1-01-FT-R-29	20
1-01-FT-L-30	20
1-01-FT-R-30	20
1-01-FT-01-21	20,87
1-01-FT-01-22	20,87
1-01-FT-01-23	20,87
1-01-FT-01-24	20,87
1-01-FT-01-25	20,87
1-01-FT-01-26	20,87
1-01-FT-01-27	20,87
1-01-FT-01-28	20,87
1-01-FTT-L-18	19
1-01-FTT-R-18	19
1-01-FTT-L-19	19
1-01-FTT-R-19	19

Product Code Feet & Adaptors	Page no.
1-01-FTT-L-20	19
1-01-FTT-R-20	19
1-01-FTT-L-21	19
1-01-FTT-R-21	19
1-01-FTT-L-22	19
1-01-FTT-R-22	19
1-01-FTT-L-23	19
1-01-FTT-R-23	19
1-01-FTT-L-24	19
1-01-FTT-R-24	19
1-01-FTT-L-25	19
1-01-FTT-R-25	19
1-01-FTT-L-26	19
1-01-FTT-R-26	19
1-01-FTT-L-27	19
1-01-FTT-R-27	19
1-01-FTT-L-28	19
1-01-FTT-R-28	19
1-01-FTT-L-29	19
1-01-FTT-R-29	19
1-01-FTT-L-30	19
1-01-FTT-R-30	19
1-01-FTT-01	19,87
1-02-2N1A	23
1-02-2N1S	23
1-02-2N1T	23
1-02-2SS-01G	24
1-02-2SS-01B	24
1-02-2SS-01P	24
1-02-2SS-S	24
1-02-2SS-L	24
1-02-2SSW	32
1-02-2ST-S	24
1-02-2ST-L	24
1-02-A1A	26
1-02-A1A-01	26
1-02-A1S	26
1-02-A1S-01	26
1-02-A1T	26
1-02-A1T-01	26
1-02-A22	27
1-02-A22W	27
1-02-A2A	27
1-02-A2S	27
1-02-A2T	27
1-02-A2ZS	28
1-02-A2ZT	28
1-02-A3S-32	31
1-02-A3S-45	31

Index

Product Code Feet & Adaptors	Page no.
1-02-A3S-60	31
1-02-A3S-75	31
1-02-A3T-32	31
1-02-A3T-45	31
1-02-A3T-60	31
1-02-A3T-75	31
1-02-A4S-32	32
1-02-A4S-45	32
1-02-A4S-60	32
1-02-A4S-75	32
1-02-A4T-32	32
1-02-A4T-45	32
1-02-A4T-60	32
1-02-A4T-75	32
1-02-A5A	28
1-02-A5S	28
1-02-A5T	28
1-02-A5ZS	28
1-02-A5ZT	28
1-02-A6A	30
1-02-A6S	30
1-02-A6T	30
1-02-A6TOFF	30
1-02-A7S	29
1-02-A7T	29
1-02-A8S	26
1-02-A8T	26
1-02-A8ZS	27
1-02-A8ZT	27
1-02-A10S	30
1-02-A10T	30
1-02-A11ZS	29
1-02-A11ZT	29
1-02-A12S	28
1-02-A12S-01	28
1-02-A12ZS	29
1-02-A12ZT	29
1-02-A13T	30
1-02-A14S	31
1-02-A14T	31
1-02-A15S	31
1-02-A15T	31
1-02-B200	25
1-02-B400	25
1-02-BO200	25
1-02-BO250	25
1-02-BO400	25
1-02-BT200	25
1-02-BT400	25

Product Code Feet & Adaptors	Page no.
1-02-C300	26
1-02-C355	26
1-02-DB1	24
1-02-T2A1S-210	24
1-02-T2A1S-410	24
1-02-T2A1T-210	24
1-02-T2A1T-410	24
1-02-T3S-45	24
1-02-T3S-60	24
1-02-T3S-75	24
1-02-T3T-45	24
1-02-T3T-60	24
1-02-T3T-75	24
1-02-TA1	28
1-02-TO1	26
1-02-TO2	26

Product Code Kneejoints & Hip Joints	Page no.
2-01-3S40	47
2-01-3S40-01	47
2-01-4S17	46
2-01-4S17-01	46
2-01-4S1S	48
2-01-4S1S-01	48
2-01-4S1S-02	48
2-01-4S1T	48
2-01-4S1T-01	48
2-01-4S1T-02	48
2-01-4S21	45
2-01-4S21-01	45
2-01-4S22	46
2-01-4S23	45
2-01-4S23-01	45
2-01-4S23-02	45
2-01-4S2S	48
2-01-4S2S-01	48
2-01-4S2S-02	48
2-01-4S2T	48
2-01-4S2T-01	48
2-01-4S2T-02	48
2-01-4S33	46
2-01-4S33-01	46
2-01-4S34	46
2-01-4S3S	44
2-01-4S3T	44
2-01-4S40	46
2-01-4S40-01	46
2-01-4S40J	47
2-01-4S40J-01	47
2-01-A41	43
2-01-A41=LA	43
2-01-A41=SP	43
2-01-A41-01	43
2-01-A41=LA-01	43
2-01-A41=SP-01	43
2-01-KA5	44
2-01-S400	42
2-01-S400-01	42
2-01-S400=LA	42
2-01-S400=LA-01	42
2-01-S400=SP	42
2-01-S400=SP-01	42
2-01-S400L	41
2-01-S400L-01	41
2-01-S400L-02	41
2-01-S400L=LA	41
2-01-S400L=LA-01	41

Index

Product Code Kneejoints & Hip Joints	Page no.
2-01-S400L=SP	41
2-01-S400L=SP-01	41
2-01-S500	40
2-01-S500-01	40
2-01-S500=LA	40
2-01-S500=LA-01	40
2-01-S500=SP	40
2-01-S500=SP-01	40
2-01-S600	39
2-01-S600-01	39
2-01-S600=LA	39
2-01-S600=LA-01	39
2-01-S600=SP	39
2-01-S600=SP-01	39
2-01-S700	38
2-01-S700-01	38
2-01-S700=LA	38
2-01-S700=LA-01	38
2-01-S700=SP	38
2-01-S700=SP-01	38
2-01-S700L	37
2-01-S700L-01	37
2-01-S700L-02	37
2-01-S700L=LA	37
2-01-S700L=LA-01	37
2-01-S700L=SP	37
2-01-S700L=SP-01	37
2-02-H1	49
2-02-H1-01	49
2-02-H1S	49
2-02-H1S-01	49
2-02-H1SL-L	50
2-02-H1SL-R	50
2-02-H1SL-01	50

Product Code Suspension Systems	Page no.
3-01-3010C-18	57
3-01-3010C-20	57
3-01-3010C-22	57
3-01-3010C-23.5	57
3-01-3010C-25	57
3-01-3010C-26.5	57
3-01-3010C-28	57
3-01-3010C-30	57
3-01-3010C-32	57
3-01-3010C-34	57
3-01-3010C-36	57
3-01-3010L-18	57
3-01-3010L-20	57
3-01-3010L-22	57
3-01-3010L-23.5	57
3-01-3010L-25	57
3-01-3010L-26.5	57
3-01-3010L-28	57
3-01-3010L-30	57
3-01-3010L-32	57
3-01-3010L-34	57
3-01-3010L-36	57
3-01-3020C-18	58
3-01-3020C-20	58
3-01-3020C-22	58
3-01-3020C-23.5	58
3-01-3020C-25	58
3-01-3020C-26.5	58
3-01-3020C-28	58
3-01-3020C-30	58
3-01-3020C-32	58
3-01-3020C-34	58
3-01-3020C-36	58
3-01-3020L-18	58
3-01-3020L-20	58
3-01-3020L-22	58
3-01-3020L-23.5	58
3-01-3020L-25	58
3-01-3020L-26.5	58
3-01-3020L-28	58
3-01-3020L-30	58
3-01-3020L-32	58
3-01-3020L-34	58
3-01-3020L-36	58
3-01-3030C-21	61
3-01-3030C-22	61
3-01-3030C-23.5	61
3-01-3030C-25	61
3-01-3030C-26.5	61

Product Code Suspension Systems	Page no.
3-01-3030C-28	61
3-01-3030C-30	61
3-01-3030C-32	61
3-01-3030C-34	61
3-01-3030C-36	61
3-01-3030C-38	61
3-01-3030C-40	61
3-01-3030C-42	61
3-01-3030C-44	61
3-01-3030L-21	61
3-01-3030L-22	61
3-01-3030L-23.5	61
3-01-3030L-25	61
3-01-3030L-26.5	61
3-01-3030L-28	61
3-01-3030L-30	61
3-01-3030L-32	61
3-01-3030L-34	61
3-01-3030L-36	61
3-01-3030L-38	61
3-01-3030L-40	61
3-01-3030L-42	61
3-01-3030L-44	61
3-01-3070C-30	59
3-01-3070C-32	59
3-01-3070C-34	59
3-01-3070C-36	59
3-01-3070C-38	59
3-01-3070C-40	59
3-01-3070C-45	59
3-01-3070L-30	59
3-01-3070L-32	59
3-01-3070L-34	59
3-01-3070L-36	59
3-01-3070L-38	59
3-01-3070L-40	59
3-01-3070L-45	59
3-01-3090L-30	60
3-01-3090L-32	60
3-01-3090L-34	60
3-01-3090L-36	60
3-01-3090L-38	60
3-01-3090L-40	60
3-01-3090L-45	60
3-01-3180M	52
3-01-3180-MD	52
3-01-AKSS-L-M	66
3-01-AKSS-R-M	66
3-01-AKSS-L-L	66

Index

Product Code Suspension Systems	Page no.
3-01-AKSS-R-L	66
3-01-AKSS-L-XL	66
3-01-AKSS-R-XL	66
3-01-AKSS-L-XXL	66
3-01-AKSS-R-XXL	66
3-01-BSL-00	63
3-01-BSL-01	63
3-01-DSTC	53
3-01-ECSS-S-BK	67
3-01-ECSS-M-BK	67
3-01-ECSS-L-BK	67
3-01-ECSS-S-TN	67
3-01-ECSS-M-TN	67
3-01-ECSS-L-TN	67
3-01-EF-S	64
3-01-EF-M	64
3-01-EF-L	64
3-01-EF-XL	64
3-01-KB-S	66
3-01-KB-M	66
3-01-KB-L	66
3-01-LCL	63
3-01-MAG30-18	54
3-01-MAG30-20	54
3-01-MAG30-22	54
3-01-MAG30-23.5	54
3-01-MAG30-25	54
3-01-MAG30-26.5	54
3-01-MAG30-28	54
3-01-MAG30-30	54
3-01-MAG30-32	54
3-01-MAG30-34	54
3-01-MAG30-36	54
3-01-MAG40-18	53
3-01-MAG40-20	53
3-01-MAG40-22	53
3-01-MAG40-23.5	53
3-01-MAG40-25	53
3-01-MAG40-26.5	53
3-01-MAG40-28	53
3-01-MAG40-30	53
3-01-MAG40-32	53
3-01-MAG40-34	53
3-01-MAG40-36	53
3-01-MAG50-30	56
3-01-MAG50-32	56
3-01-MAG50-34	56
3-01-MAG50-36	56
3-01-MAG50-38	56

Product Code Suspension Systems	Page no.
3-01-MAG50-40	56
3-01-MAG50-45	56
3-01-MAG60-30	55
3-01-MAG60-32	55
3-01-MAG60-34	55
3-01-MAG60-36	55
3-01-MAG60-38	55
3-01-MAG60-40	55
3-01-MAG60-45	55
3-01-SCL	62
3-01-SCL-P	62
3-01-SSL	62
3-01-SSL-P	62
3-01-TSL	62
3-01-TSL-S	62
3-01-V1	66
3-01-V2	65
3-01-V4	65
3-01-V5	65
3-01-V6	65
3-01-V7	64
3-01-V8	64

Product Code Cosmetic Tools & Materials	Page no.
4-01-AKSTO	72
4-01-BKSTO	72
4-01-FA-L-36	70
4-01-FA-R-36	70
4-01-FA-L-40	70
4-01-FA-R-40	70
4-01-FA-L-44	70
4-01-FA-R-44	70
4-01-FB	71
4-01-FBE	72
4-01-FBP	71
4-01-FKD-L-36	70

Product Code Cosmetic Tools & Materials	Page no.
4-01-FKD-R-36	70
4-01-FKD-L-40	70
4-01-FKD-R-40	70
4-01-FKD-L-44	70
4-01-FKD-R-44	70
4-01-FSD	72
4-02-CDND	77
4-02-CM	78
4-02-CMH	78
4-02-CSM-10	75
4-02-CSM-12	75
4-02-CSM-15	75
4-02-C-TB-3-10	81
4-02-C-TB-3-50	81
4-02-C-TB-4-10	81
4-02-C-TB-4-50	81
4-02-C-TB-5-10	81
4-02-C-TB-5-50	81
4-02-C-TB-6-10	81
4-02-C-TB-6-50	81
4-02-C-TB-8-10	81
4-02-C-TB-8-50	81
4-02-C-TB-10-10	81
4-02-C-TB-10-50	81
4-02-DDS	86
4-02-DND	77
4-02-ESM-6	75
4-02-ESM-9	75
4-02-ESM-12	75
4-02-GP-BL	87
4-02-GP-RE	87
4-02-LOC242	89
4-02-LOC271	89
4-02-MAD	89
4-02-MP-2-ARMY	74
4-02-MP-2-BLWT	74
4-02-MP-2-MCOL	74
4-02-MP-2-PKBK	74
4-02-MP-2-PUBK	74
4-02-MP-2-WTBK	74
4-02-MP-2-YLBK	74
4-02-MP-4-ARMY	74
4-02-MP-4-BLWT	74
4-02-MP-4-MCOL	74
4-02-MP-4-PKBK	74
4-02-MP-4-PUBK	74
4-02-MP-4-WTBK	74
4-02-MP-4-YLBK	74
4-02-NSP-3-10	80

Index

Product Code Cosmetic Tools & Materials	Page no.
4-02-NSP-4-10	80
4-02-NSP-5-10	80
4-02-NSP-6-10	80
4-02-NSP-8-10	80
4-02-NSP-10-10	80
4-02-NSP-12-10	80
4-02-NSP-3-50	80
4-02-NSP-4-50	80
4-02-NSP-5-50	80
4-02-NSP-6-50	80
4-02-NSP-8-50	80
4-02-NSP-10-50	80
4-02-NSP-12-50	80
4-02-NSP-3K-1-30	81
4-02-NSP-3K-3-30	81
4-02-OAK	75
4-02-PAKCB	76
4-02-PARALEX-TB-3-10	82
4-02-PARALEX-TB-3-50	82
4-02-PARALEX-TB-4-10	82
4-02-PARALEX-TB-4-50	82
4-02-PARALEX-TB-5-10	82
4-02-PARALEX-TB-5-50	82
4-02-PARALEX-TB-6-10	82
4-02-PARALEX-TB-6-50	82
4-02-PARALEX-TB-8-10	82
4-02-PARALEX-TB-8-50	82
4-02-PARALEX-TB-10-10	82
4-02-PARALEX-TB-10-50	82
4-02-PBFL	92
4-02-PBHE	92
4-02-PBHR	92
4-02-PBRD	92
4-02-PB-75	79
4-02-PB-150	79
4-02-PBOLT-1	88
4-02-PBOLT-2	88
4-02-PBOLT-3	88
4-02-PBOLT-4	88
4-02-PC	88
4-02-PDND	77
4-02-PE-2	74
4-02-PE-3	74
4-02-PE-4	74
4-02-PE-5	74
4-02-PF-2	73
4-02-PF-3	73
4-02-PF-4	73
4-02-PF-5	73

Product Code Cosmetic Tools & Materials	Page no.
4-02-PF-6	73
4-02-PIG-CA-100g	78
4-02-PIG-CA-1kg	78
4-02-PIG-BL-100g	78
4-02-PIG-BL-1kg	78
4-02-PIG-BW-100g	78
4-02-PIG-BW-1kg	78
4-02-PIG-RE-100g	78
4-02-PIG-RE-1kg	78
4-02-PIG-WT-100g	78
4-02-PIG-WT-1kg	78
4-02-PIG-YL-100g	78
4-02-PIG-YL-1kg	78
4-02-PL-1	88
4-02-PL-2	88
4-02-PP-2	74
4-02-PP-3	74
4-02-PP-4	74
4-02-PP-5	74
4-02-PUFA	79
4-02-PUFB	79
4-02-PVAB-XS	85
4-02-PVAB-S	85
4-02-PVAB-M	85
4-02-PVAB-L	85
4-02-PVAB-XL	85
4-02-PVA-E50	85
4-02-PVA-E100	85
4-02-PVA-P100	85
4-02-PZ-3	73
4-02-PZ-6	73
4-02-PZ-12	73
4-02-SD01	90
4-02-SD02I	90
4-02-SD02M	90
4-02-SD03I	90
4-02-SD03M	90
4-02-SG-14.2	86
4-02-SG-85	86
4-02-SL01-040	90
4-02-SL01-060	90
4-02-SL01-080	90
4-02-SL01-100	90
4-02-SL02-040	90
4-02-SL02-060	90
4-02-SL02-080	90
4-02-SL02-100	90
4-02-SL03-060	90
4-02-SL03-080	90

Product Code Cosmetic Tools & Materials	Page no.
4-02-SL03-100	90
4-02-SP-A	86
4-02-SP-B	86
4-02-STO-COT-060	84
4-02-STO-COT-090	84
4-02-STO-COT-130	84
4-02-STO-COT-150	84
4-02-STO-COT-250	84
4-02-STO-COT-300	84
4-02-STO-FG-080	84
4-02-STO-FG-100	84
4-02-STO-FG-130	84
4-02-STO-FG-180	84
4-02-STO-FG-250	84
4-02-STO-PER-080	84
4-02-STO-PER-100	84
4-02-STO-PER-120	84
4-02-STO-PER-130	84
4-02-STO-PER-150	84
4-02-STO-PER-180	84
4-02-STO-PER-250	84
4-02-UG-1	73
4-02-UG-2	73
4-02-ULTC-TB-3-10	83
4-02-ULTC-TB-3-50	83
4-02-ULTC-TB-4-10	83
4-02-ULTC-TB-4-50	83
4-02-ULTC-TB-5-10	83
4-02-ULTC-TB-5-50	83
4-02-ULTC-TB-6-10	83
4-02-ULTC-TB-6-50	83
4-02-ULTC-TB-8-10	83
4-02-ULTC-TB-8-50	83
4-02-ULTC-TB-10-10	83
4-02-ULTC-TB-10-50	83
4-02-ULTC-TP-1-30	83
4-02-ULTC-TP-3-30	83
4-02-ULTC-TP-3-100	83
4-02-WC-IM-L	91
4-02-WC-IM-M	91
4-02-WC-IM-S	91
4-02-WC-MT-L	91
4-02-WC-MT-M	91
4-02-WC-MT-S	91
4-02-WH-40	91
4-02-WH-45	91
4-02-WH-50	91
4-02-WS-150	89
4-02-WS-160	89

Index

Product Code Orthotics & Rehabilitation	Page no.
5-01-AFO001-L-S	94
5-01-AFO001-R-S	94
5-01-AFO001-L-M	94
5-01-AFO001-R-M	94
5-01-AFO001-L-L	94
5-01-AFO001-R-L	94
5-01-AFO001-L-XL	94
5-01-AFO001-R-XL	94
5-01-AFO002-L-S	94
5-01-AFO002-R-S	94
5-01-AFO002-L-M	94
5-01-AFO002-R-M	94
5-01-AFO002-L-L	94
5-01-AFO002-R-L	94
5-01-AFO002-L-XL	94
5-01-AFO002-R-XL	94
5-01-COL001-2-S	99
5-01-COL001-2-M	99
5-01-COL001-2-L	99
5-01-COL001-2-XL	99
5-01-COL001-3-S	99
5-01-COL001-3-M	99
5-01-COL001-3-L	99
5-01-COL001-3-XL	99
5-01-COL001-4-S	99
5-01-COL001-4-M	99
5-01-COL001-4-L	99
5-01-COL001-4-XL	99
5-01-COL001-5-S	99
5-01-COL001-5-M	99
5-01-COL001-5-L	99
5-01-COL001-5-XL	99
5-01-FL-XS	101
5-01-FL-S	101
5-01-FL-M	101
5-01-FL-L	101
5-01-FL-XL	101
5-01-HCC-S	102
5-01-HCC-M	102
5-01-HCC-L	102
5-01-HJ	96
5-01-MAJ001-P	98
5-01-MAJ001-A	98
5-01-MAJ001-HD	98
5-01-MAJKIT-P	98
5-01-MAJKIT-A	98
5-01-MAJKIT-HD	98
5-01-OFFS-SS	100
5-01-PFSS-S	96

Product Code Orthotics & Rehabilitation	Page no.
5-01-PFSS-M	96
5-01-PFSS-L	96
5-01-QL-S	102
5-01-QL-M	102
5-01-QL-L	102
5-01-RB01	97
5-01-RFD100	95
5-01-RFS100	95
5-01-RFS100DM	95
5-01-RLOC-AL	100
5-01-RLOC-SS	100
5-01-RTB-BK	96
5-01-RTB-WT	96
5-01-SB01-S	97
5-01-SB01-M	97
5-01-SB01-L	97
5-01-SB01-XL	97
5-01-SLOC-AL	100
5-01-SLOC-SS	100
5-01-TS-S	102
5-01-TS-M	102
5-01-TS-L	102
5-01-V3CJ	98
5-01-V3CJ-01	98

Product Code Upperlimbs modular prosthesis	Page no.
6-01-EJ1-CA	104
6-01-EJ1-DKBW	104
6-01-EJ2-L	104
6-01-EJ2-R	104
6-01-SJ1-L	104
6-01-SJ1-R	104
6-02-MH1-IM-L-S	105
6-02-MH1-IM-R-S	105
6-02-MH1-IM-L-M	105
6-02-MH1-IM-R-M	105
6-02-MH1-IM-L-L	105
6-02-MH1-IM-R-L	105
6-02-MH1-IM-L-XL	105
6-02-MH1-IM-R-XL	105
6-02-MH1-MT-L-S	105
6-02-MH1-MT-R-S	105
6-02-MH1-MT-L-M	105
6-02-MH1-MT-R-M	105
6-02-MH1-MT-L-L	105
6-02-MH1-MT-R-L	105
6-02-MH1-MT-L-XL	105
6-02-MH1-MT-R-XL	105

Pediatric



SACH foot (pediatric)

0-01-FTT



0-01-FTT

Product Code	Side	Foot Length	Heel Height	Max. Load
0-01-FTT-L-12	Left	12cm	10mm	60kg
0-01-FTT-R-12	Right	12cm	10mm	60kg
0-01-FTT-L-13	Left	13cm	10mm	60kg
0-01-FTT-R-13	Right	13cm	10mm	60kg
0-01-FTT-L-14	Left	14cm	10mm	60kg
0-01-FTT-R-14	Right	14cm	10mm	60kg
0-01-FTT-L-15	Left	15cm	10mm	60kg
0-01-FTT-R-15	Right	15cm	10mm	60kg
0-01-FTT-L-16	Left	16cm	10mm	60kg
0-01-FTT-R-16	Right	16cm	10mm	60kg
0-01-FTT-L-17	Left	17cm	10mm	60kg
0-01-FTT-R-17	Right	17cm	10mm	60kg

Male socket adaptor (pediatric)

0-02-A1SC/0-02-A1AC

New!



0-02-A1SC

0-02-A1AC

Product Code	Weight	Material	Max. Load
0-02-A1SC	55g	SS	45kg
0-02-A1AC	33g	AL	45kg

Female socket adaptor (pediatric)

0-02-A5AC



0-02-A5AC

New!

Product Code	Weight	Material	Max. Load
0-02-A5AC	42g	AL	45kg

Female tube adaptor (pediatric)

0-02-A6AC



0-02-A6AC

New!

Product Code	Weight	Material	Max. Load
0-02-A6AC	36g	AL	45kg

Angled tube adaptor (pediatric)

0-02-A10AC



0-02-A10AC

New!

Product Code	Weight	Material	Max. Load
0-02-A10AC	50g	AL	45kg

Attachment block (pediatric)

0-02-A22C



0-02-A22C

New!

Product Code	Weight	Material	Max. Load
0-02-A22C	98g	-	45kg

SACH ankle adaptor (pediatric)

0-02-2N1SC/0-02-2N1TC

New!



0-02-2N1SC / 0-02-2N1TC

Product Code	Weight	Material	Max. Load
0-02-2N1TC	39g	TI	45kg
0-02-2N1SC	62g	SS	45kg

Tube (pediatric)

0-02-BC150/0-02-BC300

New!



0-02-BC150



0-02-BC300

Product Code	Length	Diameter	Weight	Material	Max. Load
0-02-BC150	150mm	Ø22mm	55g	AL	45kg
0-02-BC300	300mm	Ø22mm	125g	AL	45kg

Single axis kneejoint (pediatric)

0-00-4S30C

New!

0-00-4S30C

Single axis kneejoint (pediatric)

- Lightweight aluminum single axis
- Pediatric kneejoint with proximal adjustment pyramid and distal tube clamp
- Proximal connection: adjustment pyramid
- Distal connection: tube clamp Ø 22mm



0-00-4S30C

Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S30C	160g	AL	45kg	I

Pediatric 4 Bar kneejoint (AK/KD)

0-00-4S41C

New!

0-00-4S41C

Pediatric 4 Bar kneejoint (AK/KD)

- The kneejoint is light weight 4 bar kneejoint for pediatric
- Optional adaptors are pyramid adaptor (for AK)
- 3 prong adaptor (for KD)
- If you like to order maintenance part, please feel free to contact us



0-00-4S41C



0-00-4S41C



0-00-4S41C=LA

Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S41C	320g	AL	55kg	II III
0-00-4S41C=LA	320g	AL	55kg	II III

Accessories



(AK)

0-00-4S41C-01



(KD)

0-00-4S41C=LA-01

Product Code	Material	Description
0-00-4S41C-01	SS	Pyramid adaptor(AK)
0-00-4S41C=LA-01	SS	Anchor plate(KD)

Pediatric 4 Bar kneejoint with manual lock (AK/KD)

0-00-4S40C

New!



0-00-4S40C

0-00-4S40C

Pediatric 4 Bar kneejoint with manual lock (AK/KD)

- The kneejoint is light weight manual locking kneejoint for pediatric with very low functional level
- Optional adaptors are pyramid adaptor (for AK)
- 3 prong adaptor (for KD)
- If you like to order maintenance part, please feel free to contact us



0-00-4S40C




0-00-4S40C=LA

Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S40C	387g	AL	55kg	II III
0-00-4S40C=LA	387g	AL	55kg	II III

Accessories



0-00-4S40C-01

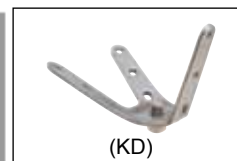

0-00-4S40C=LA-01



0-00-4S40C-02

Product Code	Material	Description
0-00-4S40C-01	SS	Pyramid adaptor(AK)
0-00-4S40C=LA-01	SS	Anchor plate(KD)
0-00-4S40C-02	-	Lock set

Pediatric aluminium polycentric mechanical hip joint

0-00-H1C



0-00-H1C

0-00-H1C Pediatric aluminium polycentric mechanical hip joint

- lightweight design

Product Code	Weight	Material	Max. Load
0-00-H1C	550g	AL	55kg



0-00-H1C-01


Product Code	Description
0-00-H1C-01	Anchor plate

Above knee foam cover (pediatric)

0-01-FAC



0-01-FAC

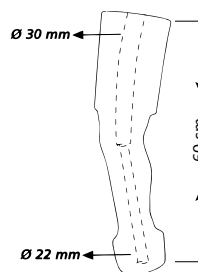

0-01-FAC Above knee foam cover (pediatric)

- Durable foam
- One size
- 30° knee flexion
- Preshaped
- Center Hole
- Flame Resistant
- Skin Color



Product Code	Calf Circum.	Hole diameter	Length
0-01-FAC	30-22cm	22-30mm	60cm


Foam covers are available in vacuum pack

Below knee foam cover (pediatric)

0-01-FBC



0-01-FBC

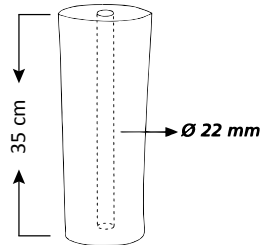

0-01-FBC

Below knee foam cover (pediatric)

- Durable foam
- One size
- Unshaped
- Center Hole
- Flame Resistant
- Skin Color

Product Code	Hole diameter	Length
0-01-FBC	22mm	35cm

Foam covers are available in vacuum pack

Ring lock system - Straight (pediatric)

0-01-RLOCC-AL



0-01-RLOCC-AL

Product Code	Length	Width	Thickness	Weight	Description
0-01-RLOCC-AL	740mm	12mm	5mm	400g	Aluminum upright bar

Ring lock system - Straight (pediatric)

0-01-RLOCC-SS



0-01-RLOCC-SS

Product Code	Length	Width	Thickness	Weight	Description
0-01-RLOCC-SS	460mm	12mm	3mm	350g	Stainless steel upright bar

Offset orthotic joint (pediatric)

0-01-OFFSC-SS



0-01-OFFSC-SS

Product Code	Length	Width	Thickness	Weight	Description
0-01-OFFSC-SS	850mm	16mm	4mm	1100g	Offset orthotic joint (stainless steel)

Feet and Adaptors



Aluminium BK kit



1-02-A1A
Male socket adaptor
page 26



1-02-A2A
Male socket adaptor
page 27



1-02-A6A
Female tube adaptor
page 30



1-02-BO250
Tube with female pyramid adaptor
page 25



1-02-2N1A
SACH ankle adaptor
page 23



1-01-FTT
SACH foot (open toe)
page 19

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Suspension
Systems

Cosmetic, Tools
& Materials

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

SACH foot (pediatric)

0-01-FTT



Product Code	Side	Foot Length	Heel Height	Max. Load
0-01-FTT-L-12	Left	12cm	10mm	60kg
0-01-FTT-R-12	Right	12cm	10mm	60kg
0-01-FTT-L-13	Left	13cm	10mm	60kg
0-01-FTT-R-13	Right	13cm	10mm	60kg
0-01-FTT-L-14	Left	14cm	10mm	60kg
0-01-FTT-R-14	Right	14cm	10mm	60kg
0-01-FTT-L-15	Left	15cm	10mm	60kg
0-01-FTT-R-15	Right	15cm	10mm	60kg
0-01-FTT-L-16	Left	16cm	10mm	60kg
0-01-FTT-R-16	Right	16cm	10mm	60kg
0-01-FTT-L-17	Left	17cm	10mm	60kg
0-01-FTT-R-17	Right	17cm	10mm	60kg

SACH foot (open toe)

1-01-FTT



Product Code	Side	Foot Length	Heel Height	Max. Load
1-01-FTT-L-18	Left	18cm	10mm	60kg
1-01-FTT-R-18	Right	18cm	10mm	60kg
1-01-FTT-L-19	Left	19cm	10mm	60kg
1-01-FTT-R-19	Right	19cm	10mm	60kg
1-01-FTT-L-20	Left	20cm	10mm	90kg
1-01-FTT-R-20	Right	20cm	10mm	90kg
1-01-FTT-L-21	Left	21cm	10mm	90kg
1-01-FTT-R-21	Right	21cm	10mm	90kg
1-01-FTT-L-22	Left	22cm	10mm	95kg
1-01-FTT-R-22	Right	22cm	10mm	95kg
1-01-FTT-L-23	Left	23cm	10mm	95kg
1-01-FTT-R-23	Right	23cm	10mm	95kg
1-01-FTT-L-24	Left	24cm	10mm	95kg
1-01-FTT-R-24	Right	24cm	10mm	95kg
1-01-FTT-L-25	Left	25cm	10mm	95kg
1-01-FTT-R-25	Right	25cm	10mm	95kg
1-01-FTT-L-26	Left	26cm	10mm	100kg
1-01-FTT-R-26	Right	26cm	10mm	100kg
1-01-FTT-L-27	Left	27cm	10mm	100kg
1-01-FTT-R-27	Right	27cm	10mm	100kg
1-01-FTT-L-28	Left	28cm	10mm	100kg
1-01-FTT-R-28	Right	28cm	10mm	100kg
1-01-FTT-L-29	Left	29cm	10mm	100kg
1-01-FTT-R-29	Right	29cm	10mm	100kg
1-01-FTT-L-30	Left	30cm	10mm	100kg
1-01-FTT-R-30	Right	30cm	10mm	100kg

1-01-FTT-01

connection plate for FTT

(one size, suitable for 1-01-FTT, size 20cm-30cm)

Product Code

1-01-FTT-01

Regal Feet and Adaptors

Single axis foot

1-01-FT



Product Code	Side	Foot Length	Heel Height	Max. Load
1-01-FT-L-21	Left	21cm	10mm	80kg
1-01-FT-R-21	Right	21cm	10mm	80kg
1-01-FT-L-22	Left	22cm	10mm	80kg
1-01-FT-R-22	Right	22cm	10mm	80kg
1-01-FT-L-23	Left	23cm	10mm	80kg
1-01-FT-R-23	Right	23cm	10mm	80kg
1-01-FT-L-24	Left	24cm	10mm	80kg
1-01-FT-R-24	Right	24cm	10mm	80kg
1-01-FT-L-25	Left	25cm	10mm	80kg
1-01-FT-R-25	Right	25cm	10mm	80kg
1-01-FT-L-26	Left	26cm	10mm	90kg
1-01-FT-R-26	Right	26cm	10mm	90kg
1-01-FT-L-27	Left	27cm	10mm	90kg
1-01-FT-R-27	Right	27cm	10mm	90kg
1-01-FT-L-28	Left	28cm	10mm	90kg
1-01-FT-R-28	Right	28cm	10mm	90kg
1-01-FT-L-29	Left	29cm	10mm	90kg
1-01-FT-R-29	Right	29cm	10mm	90kg
1-01-FT-L-30	Left	30cm	10mm	90kg
1-01-FT-R-30	Right	30cm	10mm	90kg



1-01-FT-01
connection cap for Regal FT

Product Code	Size
1-01-FT-01-21	21cm
1-01-FT-01-22	22cm
1-01-FT-01-23	23cm
1-01-FT-01-24	24cm
1-01-FT-01-25	25cm
1-01-FT-01-26	26cm
1-01-FT-01-27	27cm
1-01-FT-01-28	28cm

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

DYNA foot

1-01-DYNA

New!

1-01-DYNA

DYNA foot with fiber glass keel

- Special fiber composite and smart design to provide good energy restoring function
- An economy version compare to carbon fiber composite



1-01-DYNA

Product Code	Side	Foot length	Max. Load
1-01-DYNA-L-22	Left	22cm	75kg
1-01-DYNA-R-22	Right	22cm	75kg
1-01-DYNA-L-23	Left	23cm	75kg
1-01-DYNA-R-23	Right	23cm	75kg
1-01-DYNA-L-24	Left	24cm	75kg
1-01-DYNA-R-24	Right	24cm	75kg
1-01-DYNA-L-25	Left	25cm	80kg
1-01-DYNA-R-25	Right	25cm	80kg
1-01-DYNA-L-26	Left	26cm	80kg
1-01-DYNA-R-26	Right	26cm	80kg
1-01-DYNA-L-27	Left	27cm	80kg
1-01-DYNA-R-27	Right	27cm	80kg
1-01-DYNA-L-28	Left	28cm	80kg
1-01-DYNA-R-28	Right	28cm	80kg
1-01-DYNA-L-29	Left	29cm	80kg
1-01-DYNA-R-29	Right	29cm	80kg



1-01-DYNA-01

Product Code	Description
1-01-DYNA-01	Fixed adaptor for DYNA foot

Regal Feet and Adaptors

Pediatric

Male socket adaptor (pediatric)

0-02-A1SC/0-02-A1AC



0-02-A1SC

0-02-A1AC

New!

Product Code	Weight	Material	Max. Load
0-02-A1SC	55g	SS	45kg
0-02-A1AC	33g	AL	45kg

Female socket adaptor (pediatric)

0-02-A5AC



0-02-A5AC

New!

Product Code	Weight	Material	Max. Load
0-02-A5AC	42g	AL	45kg

Female tube adaptor (pediatric)

0-02-A6AC



0-02-A6AC

New!

Product Code	Weight	Material	Max. Load
0-02-A6AC	36g	AL	45kg

Angled tube adaptor (pediatric)

0-02-A10AC



0-02-A10AC

New!

Product Code	Weight	Material	Max. Load
0-02-A10AC	50g	AL	45kg

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Attachment block (pediatric)

0-02-A22C



0-02-A22C

New!

Product Code	Weight	Material	Max. Load
0-02-A22C	98g	-	45kg

Pediatric

Feet & Adaptors

SACH ankle adaptor

1-02-2N1T/1-02-2N1S/1-02-2N1A



1-02-2N1T

1-02-2N1S

1-02-2N1A

Product Code	Weight	Material	Max. Load
1-02-2N1T	92g	TI	125kg
1-02-2N1S	118g	SS	125kg
1-02-2N1A	98g	AL	100kg

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

SACH ankle adaptor (pediatric)

0-02-2N1SC/0-02-2N1TC



0-02-2N1SC / 0-02-2N1TC

New!

Product Code	Weight	Material	Max. Load
0-02-2N1TC	39g	TI	45kg
0-02-2N1SC	62g	SS	45kg

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Regal Feet and Adaptors

Single-axis ankle adaptor

1-02-2ST/1-02-2SS/1-02-DB1



1-02-2ST



1-02-2SS



1-02-DB1
(bearing axis)

Product Code	Size	Weight	Material	Max.Load
1-02-2ST-S	22-25cm	196g	TI	125kg
1-02-2ST-L	26-30cm	208g	TI	125kg
1-02-2SS-S	22-25cm	273g	SS	125kg
1-02-2SS-L	26-30cm	296g	SS	125kg
1-02-DB1	-	173g	AL	100kg

Heel cushion for
Single-axis ankle adaptors



1-02-2SS-01G



1-02-2SS-01B



1-02-2SS-01P

Product Code	Description	Hardness	Color
1-02-2SS-01G	Heel cushion	Hard	Green
1-02-2SS-01B	Heel cushion	Med	Black
1-02-2SS-01P	Heel cushion	Soft	Pink

Tube adaptor

1-02-T3T/1-02-T3S



75mm

60mm

45mm

1-02-T3T/
1-02-T3S

Product Code	Height	Weight	Material	Max. Load
1-02-T3T-45	45mm	48g	TI	100kg
1-02-T3T-60	60mm	56g	TI	100kg
1-02-T3T-75	75mm	67g	TI	100kg
1-02-T3S-45	45mm	84g	SS	100kg
1-02-T3S-60	60mm	99g	SS	100kg
1-02-T3S-75	75mm	120g	SS	100kg

Tube with female pyramid adaptor

1-02-T2A1T/1-02-T2A1S



1-02-T2A1T-210/
1-02-T2A1S-210



1-02-T2A1T-410/
1-02-T2A1S-410

Product Code	Length	Diameter	Weight	Tube material	Head material	Max. Load
1-02-T2A1T-210	210mm	Ø30mm	158g	AL	TI	100kg
1-02-T2A1S-210	210mm	Ø30mm	198g	AL	SS	100kg
1-02-T2A1T-410	410mm	Ø30mm	260g	AL	TI	100kg
1-02-T2A1S-410	410mm	Ø30mm	308g	AL	SS	100kg

Tube with female pyramid adaptor

1-02-BO200/1-02-BO250/1-02-BO400



1-02-BO200



1-02-BO400



1-02-BO250

Product Code	Length	Diameter	Weight	Material	Max. Load
1-02-BO200	200mm	Ø30mm	150g	AL	100kg
1-02-BO250	250mm	Ø30mm	179g	AL	100kg
1-02-BO400	400mm	Ø30mm	336g	AL	100kg

Tube with SACH ankle

1-02-BT200/1-02-BT400



1-02-BT200



1-02-BT400

Product Code	Length	Diameter	Weight	Material	Max. Load
1-02-BT200	200mm	Ø30mm	223g	AL	100kg
1-02-BT400	400mm	Ø30mm	336g	AL	100kg

Tube (pediatric)

0-02-BC150/0-02-BC300



0-02-BC150

New!



0-02-BC300

Product Code	Length	Diameter	Weight	Material	Max. Load
0-02-BC150	150mm	Ø22mm	55g	AL	45kg
0-02-BC300	300mm	Ø22mm	125g	AL	45kg

Tube

1-02-B200/1-02-B400



1-02-B200



1-02-B400

Product Code	Length	Diameter	Weight	Material	Max. Load
1-02-B200	200mm	Ø30mm	118g	AL	100kg
1-02-B400	400mm	Ø30mm	298g	AL	100kg

Regal Feet and Adaptors

Tube

1-02-C300/1-02-C355



1-02-C300



1-02-C355

Product Code	Length	Diameter	Weight	Material	Max. Load
1-02-C300	300mm	Ø30mm	115g	CF	100kg
1-02-C355	355mm	Ø30mm	133g	CF	100kg

Male socket adaptor

1-02-A1T/1-02-A1S/1-02-A1A



1-02-A1T



1-02-A1S



1-02-A1A

Product Code	Weight	Material	Max. Load
1-02-A1T	53g	TI	125kg
1-02-A1S	120g	SS	125kg
1-02-A1A	59g	AL	100kg

1-02-A1T-01/
1-02-A1S-01/
1-02-A1A-01



Lamination dummy

Product Code	Description
1-02-A1T-01	Lamination dummy
1-02-A1S-01	Lamination dummy
1-02-A1A-01	Lamination dummy

Socket adaptor

1-02-TO1/1-02-TO2



1-02-TO1



1-02-TO2

Product Code	Weight	Material	Max. Load
1-02-TO1	160g	AL	100kg
1-02-TO2	170g	AL	100kg

Female socket adaptor

1-02-A8T/1-02-A8S



1-02-A8T



1-02-A8S

Product Code	Weight	Material	Max. Load
1-02-A8T	57g	TI	100kg
1-02-A8S	96g	SS	100kg

Female socket adaptor (Rotatable)

1-02-A8ZT/1-02-A8ZS



1-02-A8ZT



1-02-A8ZS

Product Code	Weight	4 prongs Material	Head Material	Max. Load
1-02-A8ZT	65g	SS	TI	100kg
1-02-A8ZS	117g	SS	SS	100kg

PU Foam block

1-02-A22



1-02-A22



Product Code	Weight	Material	Max. Load
1-02-A22	-	Foam	-

Wooden attachment block

1-02-A22W



1-02-A22W

Product Code	Weight	Material	Max. Load
1-02-A22W	-	Wood	-

Male socket adaptor

1-02-A2T/1-02-A2S/1-02-A2A



1-02-A2T



1-02-A2S



1-02-A2A

Product Code	Weight	Material	Max. Load
1-02-A2T	43g	TI	125kg
1-02-A2S	111g	SS	125kg
1-02-A2A	56g	AL	100kg

Regal Feet and Adaptors

Pediatric

Male socket adaptor (Rotatable)

1-02-A2ZT/1-02-A2ZS



1-02-A2ZT



1-02-A2ZS

Product Code	Weight	Material	Max. Load
1-02-A2ZT	76g	TI	125kg
1-02-A2ZS	180g	SS	125kg

Female socket adaptor

1-02-A5T/1-02-A5S/1-02-A5A



1-02-A5T



1-02-A5S



1-02-A5A

Product Code	Weight	Material	Max. Load
1-02-A5T	46g	TI	125kg
1-02-A5S	76g	SS	125kg
1-02-A5A	44g	AL	100kg

Female socket adaptor (Rotatable)

1-02-A5ZT/1-02-A5ZS/1-02-TA1



1-02-A5ZT



1-02-A5ZS



1-02-TA1

Product Code	Weight	Material	Max. Load
1-02-A5ZT	90g	TI	125kg
1-02-A5ZS	142g	SS	125kg
1-02-TA1	102g	AL	100kg

Male socket adaptor

1-02-A12S



1-02-A12S



1-02-A12S-01

Product Code	Weight	Material	Max. Load
1-02-A12S	124g	SS	100kg

Product Code	Description
1-02-A12S-01	Lamination dummy

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Male socket adaptor (Rotatable)

1-02-A12ZT/1-02-A12ZS



1-02-A12ZT



1-02-A12ZS

Product Code	Weight	3 prongs Material	Head Material	Max. Load
1-02-A12ZT	154g	SS	Ti	125kg
1-02-A12ZS	178g	SS	SS	125kg

Female socket adaptor (Rotatable)

1-02-A11ZT/1-02-A11ZS



1-02-A11ZT



1-02-A11ZS

Product Code	Weight	3 prongs Material	Head Material	Max. Load
1-02-A11ZT	154g	SS	Ti	125kg
1-02-A11ZS	178g	SS	SS	125kg

Anchor plate (KD)

0-00-4S41C=LA-01/0-00-4S40C=LA-01



0-00-4S41C=LA-01/
0-00-4S40C=LA-01

New!

Product Code	Material	Description
0-00-4S41C=LA-01	SS	Anchor plate(KD)
0-00-4S40C=LA-01	SS	Anchor plate(KD)

Male tube adaptor

1-02-A7T/1-02-A7S



1-02-A7T



1-02-A7S

Product Code	Weight	Material	Max. Load
1-02-A7T	60g	Ti	100kg
1-02-A7S	96g	SS	100kg

Female tube adaptor

1-02-A6T/1-02-A6S/1-02-A6A



1-02-A6T



1-02-A6A



1-02-A6S

Product Code	Weight	Material	Max. Load
1-02-A6T	63g	TI	125kg
1-02-A6S	113g	SS	125kg
1-02-A6A	78g	AL	100kg

Offset Tube Clamp

1-02-A6TOFF



1-02-A6TOFF

1-02-A6TOFF Offset Tube Clamp

- The titanium 30mm offset tube clamp is used when extra alignment is needed
- It is also ideal when using a pin suspension lock, where prosthetic and pin alignment can be easily adjusted using the offset clamp
- This versatile clamp design also provides an offset position that can be rotated 360 degrees

Product Code	Weight	Material	Max. Load
1-02-A6TOFF	109g	TI	125kg

Angled tube adaptor

1-02-A10T/1-02-A10S



1-02-A10T



1-02-A10S

Product Code	Weight	Material	Max. Load
1-02-A10T	73g	TI	125kg
1-02-A10S	125g	SS	125kg

Rotation adaptor

1-02-A13T



1-02-A13T

1-02-A13T Rotation adaptor

- Set Screws providing the lower connection
- With the incorporation of the rotation adapter above the knee
- The shin may be rotated relative to the socket with the knee flexed
- For rotation, press the button then turn the knee
- Suitable for the trans-femoral amputees, hip disarticulation amputees

Product Code	Weight	Material	Max. Load
1-02-A13T	175g	TI	100kg

Sliding tube clamp

1-02-A14T/1-02-A14S



1-02-A14T



1-02-A14S

1-02-A14T/1-02-A14S Sliding tube clamp

- 30mm sliding tube clamp adaptor
- Permits additional sliding adjustment between the prosthetic socket and tube adaptor, even as the prosthetic is being worn
- The sliding adjustment allows for fine turning in any desired parallel position of the distal portion, independent of the angle of adjustment through the pyramid

Product Code	Weight	Material	Max. Load
1-02-A14T	180g	TI	100kg
1-02-A14S	313g	SS	100kg

Double adaptor with adjusting of length

1-02-A15T/1-02-A15S



1-02-A15T/1-02-A15S

1-02-A15T/1-02-A15S Double adaptor with adjusting of length

- Suitable for TT and TF amputee with long residual limb
- The height could be adjusted from 60mm to 75mm

Product Code	Weight	Height	Material	Max. Load
1-02-A15T	133g	60-75mm	TI	100kg
1-02-A15S	205g	60-75mm	SS	100kg

Double adaptor

1-02-A3T/1-02-A3S



75mm

60mm

45mm

32mm

1-02-A3T/
1-02-A3S

Product Code	Weight	Material	Max. Load
1-02-A3T-32	63g	TI	125kg
1-02-A3T-45	72g	TI	125kg
1-02-A3T-60	84g	TI	125kg
1-02-A3T-75	90g	TI	125kg
1-02-A3S-32	115g	SS	125kg
1-02-A3S-45	144g	SS	125kg
1-02-A3S-60	174g	SS	125kg
1-02-A3S-75	202g	SS	125kg

Double adaptor

1-02-A4T/1-02-A4S



Product Code	Weight	Material	Max. Load
1-02-A4T-32	70g	TI	125kg
1-02-A4T-45	82g	TI	125kg
1-02-A4T-60	93g	TI	125kg
1-02-A4T-75	103g	TI	125kg
1-02-A4S-32	121g	SS	125kg
1-02-A4S-45	163g	SS	125kg
1-02-A4S-60	192g	SS	125kg
1-02-A4S-75	222g	SS	125kg

Single axis ankle joint (for wooden single axis ankle)

1-02-2SSW



1-02-2SSW

New!

Product Code

1-02-2SSW

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Kneejoints and Hip Joints



Pediatric 4 Bar kneejoint (AK/KD)

0-00-4S41C

New!

0-00-4S41C
Pediatric 4 Bar kneejoint (AK/KD)

- The kneejoint is light weight 4 bar kneejoint for pediatric
- Optional adaptors are pyramid adaptor (for AK)
- 3 prong adaptor (for KD)
- If you like to order maintenance part, please feel free to contact us



0-00-4S41C



0-00-4S41C



0-00-4S41C=LA


Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S41C	320g	AL	55kg	II III
0-00-4S41C=LA	320g	AL	55kg	II III

Accessories



(AK)

0-00-4S41C-01

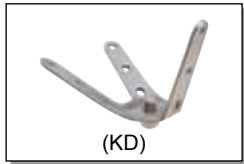

(KD)

0-00-4S41C=LA-01

Product Code	Material	Description
0-00-4S41C-01	SS	Pyramid adaptor(AK)
0-00-4S41C=LA-01	SS	Anchor plate(KD)

Pediatric 4 Bar kneejoint with manual lock (AK/KD)

0-00-4S40C

New!



0-00-4S40C

0-00-4S40C

Pediatric 4 Bar kneejoint with manual lock (AK/KD)

- The kneejoint is light weight manual locking kneejoint for pediatric with very low functional level
- Optional adaptors are pyramid adaptor (for AK)
- 3 prong adaptor (for KD)
- If you like to order maintenance part, please feel free to contact us



0-00-4S40C



0-00-4S40C=LA


Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S40C	387g	AL	55kg	II III
0-00-4S40C=LA	387g	AL	55kg	II III

Accessories



(AK)

0-00-4S40C-01

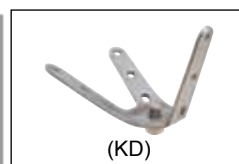

(KD)

0-00-4S40C=LA-01



Lock set

0-00-4S40C-02

Product Code	Material	Description
0-00-4S40C-01	SS	Pyramid adaptor(AK)
0-00-4S40C=LA-01	SS	Anchor plate(KD)
0-00-4S40C-02	-	Lock set

Pediatric aluminium polycentric mechanical hip joint

0-00-H1C



0-00-H1C

0-00-H1C Pediatric aluminium polycentric mechanical hip joint

- lightweight design

Product Code	Weight	Material	Max. Load
0-00-H1C	550g	AL	55kg



0-00-H1C-01

Product Code	Description
0-00-H1C-01	Anchor plate

Single axis kneejoint (pediatric)

0-00-4S30C



0-00-4S30C

New!

0-00-4S30C

Single axis kneejoint (pediatric)

- Lightweight aluminum single axis
- Pediatric kneejoint with proximal adjustment pyramid and distal tube clamp
- Proximal connection: adjustment pyramid
- Distal connection: tube clamp Ø 22mm

Product Code	Weight	Material	Max. Load	Mobility Level
0-00-4S30C	160g	AL	45kg	I

Polycentric manual locking K2 kneejoint w/ stance flexion (AK/KD) 2-01-S700L



2-01-S700L

2-01-S700L


Polycentric manual locking K2 kneejoint with stance flexion (AK/KD)

The new standard in safety knees

This NEW knee is not the friction braking knee of old. This geometrically aligned system never fails in stability

The S700L is a polycentric manual locking knee with stance flexion. The stance flexion system forces the geometry into a very stable position at heel contact until the patient rolls over the toe to release the knee into swing. The knee can be flexed in stance while providing stance stability. Elegantly simplistic design for stability. With correct alignment the knee is always in a stable position for patient confidence.

- Adjustable stance flexion from extreme stability with high stance flexion of up to 20 degrees
- Rapid response with low stance flexion
- Adjustable constant friction for swing phase resistance.
- Adjustable extension assist for swing phase
- Lock is manually released by pulling on the lanyard which is attached to the knee lock mechanism
- Lock mechanism can also be removed at a later time so that the knee can be used as a free swing knee
- Usable for either above knee or knee disarticulation
- Best suited to be used for K2
- The frame construction is made of superlight aluminum alloy
- The proximal attachments has A-P slide and Rotation adjustment
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)
- Lanyard is included with the purchase of the S700L

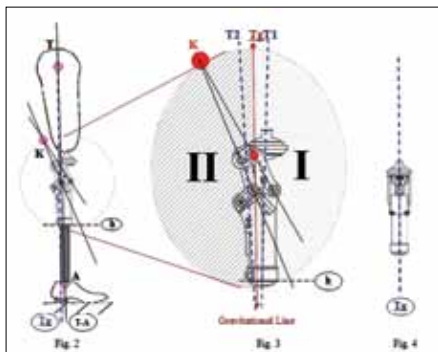

2-01-S700L

2-01-S700L=LA

2-01-S700L=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S700L	852g	AL	125kg	II III
2-01-S700L=LA	869g	AL	125kg	II III
2-01-S700L=SP	842g	AL	125kg	II III

Accessories




2-01-S700L-01



2-01-S700L=LA-01



2-01-S700L=SP-01

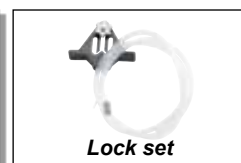

2-01-S700L-02

Product Code	Material	Description
2-01-S700L-01	AL	Pyramid adaptor(AK)
2-01-S700L=LA-01	AL	Lotus adaptor(KD)
2-01-S700L=SP-01	AL	Socket attachment plate(KD)
2-01-S700L-02	-	Lock set

Polycentric K2 kneejoint with stance flexion system (AK/KD) 2-01-S700



2-01-S700

2-01-S700


Polycentric K2 kneejoint with stance flexion system (AK/KD)

The new standard in safety knees

This NEW knee is not the friction braking knee of old. This geometrically aligned system never fails in stability

The S700 polycentric K2 knee with stance flexion system forces the geometry into a very stable position at heel contact until the patient rolls over the toe to release the knee into swing. The knee can be flexed in stance up to 20 degrees while providing stance stability. Elegantly simplistic design for stability. With correct alignment the knee is always in a stable position for patient confidence.

- Adjustable stance flexion from extreme stability with high stance flexion of up to 20 degrees to a very rapid response with low stance flexion
- Adjustable constant friction for swing phase resistance
- Adjustable extension assist for swing phase
- Usable for either above knee or knee disarticulation
- Best suited to be used for K2
- The frame construction is made of superlight aluminum alloy
- The proximal attachments has A-P slide and Rotation adjustment
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)

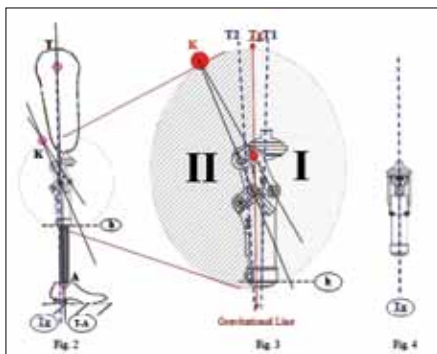

2-01-S700



2-01-S700=LA



2-01-S700=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S700	782g	AL	125kg	II III
2-01-S700=LA	799g	AL	125kg	II III
2-01-S700=SP	772g	AL	125kg	II III

Accessories



Pyramid adaptor(AK)

2-01-S700-01



Lotus adaptor(KD)

2-01-S700=LA-01



Socket attachment plate (KD)

2-01-S700=SP-01

Product Code	Material	Description
2-01-S700-01	AL	Pyramid adaptor(AK)
2-01-S700=LA-01	AL	Lotus adaptor(KD)
2-01-S700=SP-01	AL	Socket attachment plate(KD)

5 Bar polycentric pneumatic kneejoint w/ stance flexion system (AK/KD) 2-01-S600



2-01-S600


2-01-S600

5 Bar polycentric pneumatic kneejoint with stance flexion system (AK/KD)

The S600, is designed to be a very stable 5 Bar Knee with a large stance flexion angle (up to 12°). Weighing in at less than 1KG and having extensive alignment adjusting abilities give the S600 knee a smooth and natural gait pattern.

The frame construction is made of superlight aluminum alloy. All axes have built-in bearing which allow for ultra smooth walking movement. The proximal attachment has A-P slide and Rotation adjustment.

- Polycentric knee with stance flexion control for maximum safety during stance and is best suited to be used for low to high K3 level
- Practitioner adjustable from extreme stability to high stance flexion
- The S600 offers a pneumatic cylinder with independent flexion extension adjustment for swing phase
- Can be used for either above knee or knee disarticulation
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)

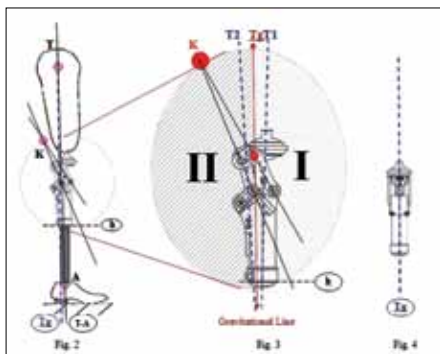

2-01-S600



2-01-S600=LA



2-01-S600=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S600	982g	AL	125kg	II III
2-01-S600=LA	999g	AL	125kg	II III
2-01-S600=SP	972g	AL	125kg	II III

Accessories



Pyramid adaptor(AK)

2-01-S600-01



Lotus adaptor(KD)

2-01-S600=LA-01



Socket attachment plate (KD)

2-01-S600=SP-01

Product Code	Material	Description
2-01-S600-01	AL	Pyramid adaptor(AK)
2-01-S600=LA-01	AL	Lotus adaptor(KD)
2-01-S600=SP-01	AL	Socket attachment plate(KD)

Polycentric 4 Bar kneejoint (AK/KD)

2-01-S500



2-01-S500

2-01-S500

Polycentric 4 Bar kneejoint (AK/KD)

The S500 knee joint is a polycentric knee with pneumatic swing control that is best suited to be used for moderate to K3 level activity and high need for safety. The frame construction is made of superlight aluminum alloy.

- Polycentric knee joint with pneumatic swing control
- Independent flexion / extension adjustment
- Built-in bearings which allow for ultra smooth walking movement
- Frame construction made of superlight aluminum alloy
- Proximal attachment has A-P slide and rotating adjustment
- Adjustable Extension assist
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)



2-01-S500



2-01-S500=LA



2-01-S500=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S500	882g	AL	125kg	II III
2-01-S500=LA	899g	AL	125kg	II III
2-01-S500=SP	872g	AL	125kg	II III

Accessories



Pyramid adaptor(AK)

2-01-S500-01



Lotus adaptor(KD)

2-01-S500=LA-01



Socket attachment plate (KD)

2-01-S500=SP-01

Product Code	Material	Description
2-01-S500-01	AL	Pyramid adaptor(AK)
2-01-S500=LA-01	AL	Lotus adaptor(KD)
2-01-S500=SP-01	AL	Socket attachment plate(KD)

Polycentric 4 Bar Linkage Kneejoint w/ Manual Lock (AK/KD)

2-01-S400L



2-01-S400L

2-01-S400L

Polycentric 4 Bar Linkage Kneejoint with Manual Lock (AK/KD)

The S400L knee joint is a manual locking, polycentric 4 bar linkage knee with friction swing phase control. The lock is manually released by pulling on the lanyard which is attached to the knee lock mechanism. The knee lock mechanism can also be removed at a later time so that the knee can be used as a free swing knee. Swing phase resistance is controlled by amount of friction applied on the axis of the knee.

- Manual locking polycentric knee
- Lanyard lock mechanism can be removed at a later time so that the knee can be used as a free swing knee
- Low activity
- High need for security
- Frame construction is made of superlight aluminum alloy
- Proximal attachment has A-P slide and rotation adjustment
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)
- Lanyard is included with the purchase of the S400L



2-01-S400L

2-01-S400L=LA

2-01-S400L=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S400L	678g	AL	125kg	II III
2-01-S400L=LA	695g	AL	125kg	II III
2-01-S400L=SP	668g	AL	125kg	II III

Accessories



Pyramid adaptor(AK)

2-01-S400L-01



Lotus adaptor(KD)

2-01-S400L=LA-01



Socket attachment plate (KD)

2-01-S400L=SP-01



Lock set

2-01-S400L-02

Product Code	Material	Description
2-01-S400L-01	AL	Pyramid adaptor(AK)
2-01-S400L=LA-01	AL	Lotus adaptor(KD)
2-01-S400L=SP-01	AL	Socket attachment plate(KD)
2-01-S400L-02	-	Lock set

Polycentric 4 Bar kneejoint (AK/KD)

2-01-S400



2-01-S400
Polycentric 4 Bar kneejoint (AK/KD)

The S400 knee joint is a polycentric knee that is best suited to be used for moderate to K3 level activity and high need for safety.

- Polycentric knee
- Constant friction swing phase control
- Frame construction is made of superlight aluminum alloy
- Proximal attachment has A-P slide and rotation adjustment
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)

2-01-S400



2-01-S400



2-01-S400=LA



2-01-S400=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-S400	650g	AL	125kg	II III
2-01-S400=LA	667g	AL	125kg	II III
2-01-S400=SP	640g	AL	125kg	II III

Accessories




Pyramid adaptor(AK)

2-01-S400-01



Lotus adaptor(KD)

2-01-S400=LA-01

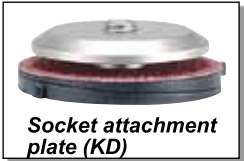

Socket attachment plate (KD)

2-01-S400=SP-01

Product Code	Material	Description
2-01-S400-01	AL	Pyramid adaptor(AK)
2-01-S400=LA-01	AL	Lotus adaptor(KD)
2-01-S400=SP-01	AL	Socket attachment plate(KD)

Polycentric 4 Bar Linkage Kneejoint (AK/KD)

2-01-A41



2-01-A41

2-01-A41

Polycentric 4 Bar Linkage Kneejoint (AK)

- Polycentric 4 bar knee
- Constant friction swing phase control
- Best suited for low activity level
- High need for security
- Frame construction is made of superlight aluminum alloy
- Proximal attachment has A-P slide and rotation adjustment
- Optional adapters: Pyramid(AK), Lotus(KD) and Socket attachment plate(KD)



2-01-A41



2-01-A41=LA



2-01-A41=SP

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-A41	559g	AL	100kg	II III
2-01-A41=LA	576g	AL	100kg	II III
2-01-A41=SP	549g	AL	100kg	II III

Accessories



Pyramid adaptor(AK)

2-01-A41-01



Lotus adaptor(KD)

2-01-A41=LA-01



Socket attachment plate (KD)

2-01-A41=SP-01

Product Code	Material	Description
2-01-A41-01	AL	Pyramid adaptor(AK)
2-01-A41=LA-01	AL	Lotus adaptor(KD)
2-01-A41=SP-01	AL	Socket attachment plate(KD)

Pneumatic 4 Bar kneejoint with integrated pyramid

2-01-KA5



2-01-KA5

New!

2-01-KA5

Pneumatic 4 Bar kneejoint with integrated pyramid

- The kneejoint is a polycentric kneejoint best suited for use with low to constant walking speed while providing a higher degree of safety
- This frame construction is made of superlight aluminum alloy and incorporates swing phase control with independent flexion/extension adjustment

Product Code	Weight	Material	Max. Load
2-01-KA5	731g	AL	100kg

4 Bar mechanical kneejoint (AK)

2-01-4S3T/2-01-4S3S



2-01-4S3T/
2-01-4S3S

2-01-4S3T/2-01-4S3S

4 Bar mechanical kneejoint (AK)

- Viscous-elastic swing
- Max.load: 100kg
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S3T	453g	TI	100kg
2-01-4S3S	662g	SS	100kg

4 Bar mechanical kneejoint (KD)

2-01-4S21



2-01-4S21

2-01-4S21

4 Bar mechanical kneejoint (KD)

- Basic level mechanical 4 linkages
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S21	870g	SS	100kg



2-01-4S21-01

Product Code	Description
2-01-4S21-01	Anchor plate

4 Bar mechanical kneejoint with lock (KD)

2-01-4S23



2-01-4S23

2-01-4S23

4 Bar mechanical kneejoint with lock (KD)

- Basic level mechanical 4 linkages
- With cable manual lock
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S23	896g	SS	100kg



2-01-4S23-02



2-01-4S23-01

Product Code	Description
2-01-4S23-01	Anchor plate
2-01-4S23-02	Lock set

Single axis lock kneejoint (AK)

2-01-4S40



2-01-4S40

Single axis lock kneejoint (AK)

- CNC machined
- Good lock design
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-4S40	278g	AL	100kg	I

Product Code	Description
2-01-4S40-01	Lock set

Single axis kneejoint with manual lock

2-01-4S33/2-01-4S17



2-01-4S33 / 2-01-4S17

Single axis kneejoint with manual lock

- with manual lock
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S33	408g	TI	125kg
2-01-4S17	586g	SS	125kg

Product Code	Description
2-01-4S33-01	Lock set
2-01-4S17-01	Lock set

Single axis kneejoint

2-01-4S34/2-01-4S22



2-01-4S34 / 2-01-4S22

Single axis kneejoint

- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S34	408g	TI	100kg
2-01-4S22	590g	SS	125kg

Single axis lock kneejoint

2-01-3S40



2-01-3S40



2-01-3S40-01

New!

2-01-3S40 Single axis lock kneejoint

- CNC machined
- Light aluminum alloy
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-3S40	394g	AL	100kg	I

Product Code	Description
2-01-3S40-01	Lock set

Single axis lock kneejoint (AK)

2-01-4S40J



2-01-4S40J




2-01-4S40J-01

New!

2-01-4S40J Single axis lock kneejoint (AK)

- CNC machined
- Good lock design
- Rotatable pyramid adaptor
- If you like to order maintenance part, please feel free to contact us

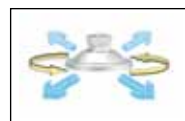

360 degree rotatable
horizontal shiftable upper pyramid

Product Code	Weight	Material	Max. Load	Mobility Level
2-01-4S40J	358g	AL	100kg	I

Product Code	Description
2-01-4S40J-01	Lock set

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Single axis kneejoint (with stance control and constant friction) 2-01-4S1T/2-01-4S1S



2-01-4S1T/
2-01-4S1S

2-01-4S1T / 2-01-4S1S Single axis kneejoint with stance control and constant friction

- with stance control and constant friction
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S1T	358g	TI	100kg
2-01-4S1S	528g	SS	100kg



2-01-4S1T-01/
2-01-4S1S-01



2-01-4S1T-02/
2-01-4S1S-02

Product Code	Description
2-01-4S1T-01	Extension Assist
2-01-4S1S-01	Extension Assist
2-01-4S1T-02	Cable
2-01-4S1S-02	Cable

Single axis kneejoint 2-01-4S2T/2-01-4S2S



4S2T/
4S2S

2-01-4S2T / 2-01-4S2S Single axis kneejoint

- with stance control and constant friction
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-01-4S2T	274g	TI	100kg
2-01-4S2S	380g	SS	100kg



2-01-4S2T-01/
2-01-4S2S-01



2-01-4S2T-02/
2-01-4S2S-02

Product Code	Description
2-01-4S2T-01	Extension Assist
2-01-4S2S-01	Extension Assist
2-01-4S2T-02	Wire
2-01-4S2S-02	Wire

Aluminium polycentric mechanical hip joint

2-02-H1



2-02-H1

2-02-H1

Aluminium polycentric mechanical hip joint

- Super lightweight design
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-02-H1	688g	AL	100kg




2-02-H1-01

Product Code	Description
2-02-H1-01	Anchor plate

Stainless steel single axis mechanical hip joint

2-02-H1S

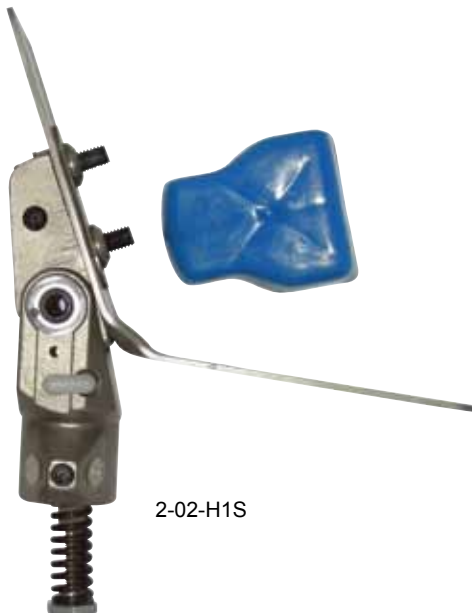

2-02-H1S

2-02-H1S

Stainless steel single axis mechanical hip joint

- Basic level mechanical single axis hip joint
- If you like to order maintenance part, please feel free to contact us

Product Code	Weight	Material	Max. Load
2-02-H1S	850g	SS	100kg



2-02-H1S-01

Product Code	Description
2-02-H1S-01	Anchor plate

Single axis mechanical hip joint with lock

2-02-H1SL

New!

2-02-H1SL
Single axis mechanical hip joint with lock
(stainelss steel & aluminium)

- For hip disarticulation and hemi-pelvectomy amputees with lower functional levels
- Higher stance phase stability, minimal swing phase control

Product Code	Side	Weight	Material	Max. Load
2-02-H1SL-L	Left	851g	SS+AL	100kg
2-02-H1SL-R	Right	851g	SS+AL	100kg



2-02-H1SL



2-02-H1SL-01

Product Code	Description
2-02-H1SL-01	Anchor plate

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty


Suspension Systems



Mechanical MAGLOCK

3-01-3180M

Pediatric
Feet & Adaptors
Kneejoints & Hip Joints
Suspension Systems
Cosmetic, Tools & Materials
Orthotics & Rehabilitation
Upperlimbs Modular Prosthesis
FAQ & Warranty

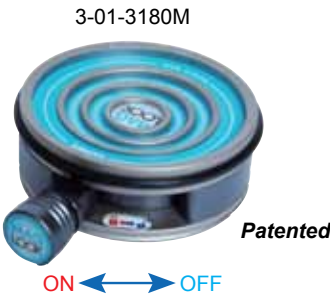

New!
3-01-3180M
Mechanical MAG Lock

- GO MAGNETIC NEW & IMPROVED MAGLOCK 2.0
- INTRODUCING THE NEWLY DESIGNED LOW PROFILE, NO BATTERY MECHANICAL
- The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3
- The MAGLOCK is developed to delivery up to 90lbs of adhesion
- The large distal target allows the patient to be 10 degrees off center and still achieve suspension
- The ultimate product for patients who have dexterity issues, issues with pin alignment, requires a caregiver and much more



Refer to "Instructional video" from our website www.regalprosthesis.com

Product Code	Weight Limit	Product Weight	Build Height	SuspensionStrength
3-01-3180M	136kg / 300lbs	278g / 0.6lbs	2.2cm / 0.87in	90lbs

Dummy for Mechanical MAGLOCK

3-01-3180-MD



New!
3-01-3180-MD
Dummy for Mechanical MAGLOCK

- Comes with different sized silicon rings for casting over a MAGLiner or when using the Distal Cap converter for a locking liner
- These rings provide the perfect spacing and alignment of the MAGLOCK
- The Fabrication Dummy comes with two silicone rings



Product Code
3-01-3180-MD

Distal Cap Converter for all pin liners - included dummy (Size: S & L) 3-01-DSTC



New!

3-01-DSTC Distal Cap Converter for all pin liners included dummy (Size: S & L)

- Cap converter with 10 mm thread for all pin liners
- The distal cap converter is required for use with other liners in order to be compatible with the MAGLOCK
- DSTC comes with two silicone rings



Product Code

3-01-DSTC

TRANSTIBIAL MAGLiner – Premier (BK) 3-01-MAG40

New!

3-01-MAG40 TRANSTIBIAL MAGLiner – Premier (BK)



Exclusively designed for the MAGLOCK

- This premium liner is designed for the health of the residual limb and ease of donning
- Manufactured with mineral oil embedded Platinum cured medical grade silicone and sheathed in a seamless durable fabric
- The matrix which is built into the liner prevents distal elongation
- 6 months warranty

Using the MAGLINER with the MAGLOCK offers the lowest profile possible



3-01-MAG40

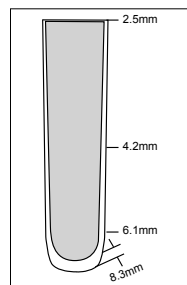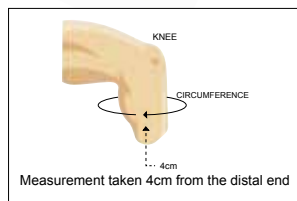

Product Code	Circumference	Length
3-01-MAG40-18	180mm	360mm
3-01-MAG40-20	200mm	360mm
3-01-MAG40-22	220mm	360mm
3-01-MAG40-23.5	235mm	360mm
3-01-MAG40-25	250mm	360mm
3-01-MAG40-26.5	265mm	360mm
3-01-MAG40-28	280mm	360mm
3-01-MAG40-30	300mm	360mm
3-01-MAG40-32	320mm	360mm
3-01-MAG40-34	340mm	360mm
3-01-MAG40-36	360mm	360mm

MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.

The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3. The MAGLOCK is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.

TRANSTIBIAL MAGLiner – Hydro (BK)

3-01-MAG30



3-01-MAG30

New!



3-01-MAG30

TRANSTIBIAL MAGLiner – Hydro (BK)

Exclusively designed for the MAGLOCK

- *Regal* Prosthesis revolutionized the liner industry with the use of Parylene® on the surface of the silicone liner. We were also the first to encapsulate the outer cloth into the liner. We have eliminated delamination. We have designed a new cloth specially for this process
- The Hydro MAGLiner's silky smooth mineral oil embedded Platinum cured medical grade silicone liner has the unique ability to repel water. This hydrophobic property makes the Hydro the easiest and fastest liner to clean and sanitize
- The Parylene® surface reduces surface friction between the liner/skin and liner/socket for better comfort and donning. The Hydro MAGLiner reduces skin shear for fragile skin and will help with blister reduction. The matrix which is built into the liner prevents distal elongation and delamination
- The Hydro Skin offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Hydro Skin liner's opaque pigment gives a natural "skin tone" to all ethnic skin types
- 6 Month warranty

Using the MAGLINER with the MAGLOCK offers the lowest profile possible

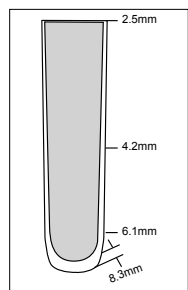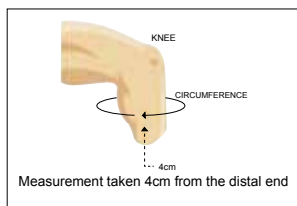

Product Code	Circumference	Length
3-01-MAG30-18	180mm	370mm
3-01-MAG30-20	200mm	370mm
3-01-MAG30-22	220mm	370mm
3-01-MAG30-23.5	235mm	370mm
3-01-MAG30-25	250mm	370mm
3-01-MAG30-26.5	265mm	370mm
3-01-MAG30-28	280mm	370mm
3-01-MAG30-30	300mm	370mm
3-01-MAG30-32	320mm	370mm
3-01-MAG30-34	340mm	370mm
3-01-MAG30-36	360mm	370mm

MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.

The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3. The MAGLOCK is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.

TRANSFEMORAL MAGLiner – Premier (AK)

3-01-MAG60



New!

3-01-MAG60



TRANSFEMORAL MAGLiner – Premier (AK)

Exclusively designed for the MAGLOCK

- This premium liner is designed for the health of the residual limb and ease of donning
- Manufactured with mineral oil embedded Platinum cured medical grade silicone and sheathed in a seamless durable fabric
- The matrix which is built into the liner prevents distal elongation
- 6 months warranty

Using the MAGLINER with the MAGLOCK offers the lowest profile possible

3-01-MAG60

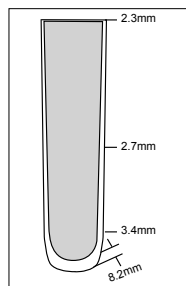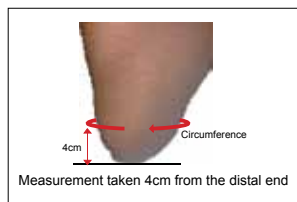

Product Code	Circumference	Length
3-01-MAG60-30	300mm	370mm
3-01-MAG60-32	320mm	370mm
3-01-MAG60-34	340mm	370mm
3-01-MAG60-36	360mm	370mm
3-01-MAG60-38	380mm	370mm
3-01-MAG60-40	400mm	370mm
3-01-MAG60-45	450mm	370mm

MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.

The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3. The MAGLOCK is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.

TRANSFEMORAL MAGLiner – Hydro (AK)

3-01-MAG50

New!



3-01-MAG50

TRANSFEMORAL MAGLiner – Hydro (AK)

Exclusively designed for the MAGLOCK

- *Regal* Prosthesis revolutionized the liner industry with the use of Parylene® on the surface of the silicone liner. We were also the first to encapsulate the outer cloth into the liner. We have eliminated delamination. We have designed a new cloth specially for this process
- The Hydro MAGLiner's silky smooth mineral oil embedded Platinum cured medical grade silicone liner has the unique ability to repel water. This hydrophobic property makes the Hydro the easiest and fastest liner to clean and sanitize
- The Parylene® surface reduces surface friction between the liner/skin and liner/socket for better comfort and donning. The Hydro MAGLiner reduces skin shear for fragile skin and will help with blister reduction. The matrix which is built into the liner prevents distal elongation and delamination
- The Hydro Skin offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Hydro Skin liner's opaque pigment gives a natural "skin tone" to all ethnic skin types
- 6 Month warranty

Using the MAGLINER with the MAGLOCK offers the lowest profile possible

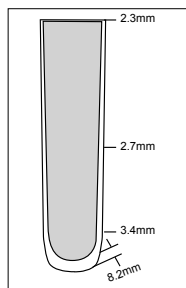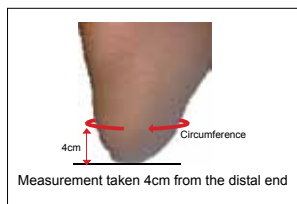

Product Code	Circumference	Length
3-01-MAG50-30	300mm	370mm
3-01-MAG50-32	320mm	370mm
3-01-MAG50-34	340mm	370mm
3-01-MAG50-36	360mm	370mm
3-01-MAG50-38	380mm	370mm
3-01-MAG50-40	400mm	370mm
3-01-MAG50-45	450mm	370mm

MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.

The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3. The MAGLOCK is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.

TRANSTIBIAL Silicone Liner – Forte (BK)

3-01-3010

3-01-3010

TRANSTIBIAL Silicone Liner (BK)

- The Forte liner is a silicone gel liner with mineral oil embedded in the silicone providing additional healing factor to the skin. The liner is made of medical grade silicone and durable seamless fabric on the outside for easy donning and doffing application. The liner comes with cushion and locking options. It is suitable for low to moderate activity patients. The matrix which is built into the liner prevents distal elongation of the liner for locking application
- 6 months warranty






Locking liner

3-01-3010L



Cushion liner

3-01-3010C

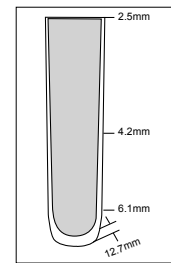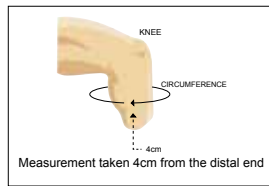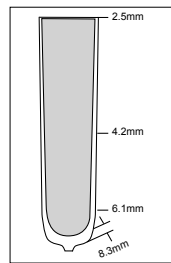

TRANSTIBIAL Locking liner

Product Code	Circumference	Length
3-01-3010L-18	180mm	370mm
3-01-3010L-20	200mm	370mm
3-01-3010L-22	220mm	370mm
3-01-3010L-23.5	235mm	370mm
3-01-3010L-25	250mm	370mm
3-01-3010L-26.5	265mm	370mm
3-01-3010L-28	280mm	370mm
3-01-3010L-30	300mm	370mm
3-01-3010L-32	320mm	370mm
3-01-3010L-34	340mm	370mm
3-01-3010L-36	360mm	370mm




TRANSTIBIAL Cushion liner

Product Code	Circumference	Length
3-01-3010C-18	180mm	360mm
3-01-3010C-20	200mm	360mm
3-01-3010C-22	220mm	360mm
3-01-3010C-23.5	235mm	360mm
3-01-3010C-25	250mm	360mm
3-01-3010C-26.5	265mm	360mm
3-01-3010C-28	280mm	360mm
3-01-3010C-30	300mm	360mm
3-01-3010C-32	320mm	360mm
3-01-3010C-34	340mm	360mm
3-01-3010C-36	360mm	360mm

3-01-3020

TRANSTIBIAL Silicone Liner (BK)

- Parasil Liner offers a superior durable design for your active patient and is available in both PIN LOCKING and CUSHION
- The silky smooth Parasil Liner is designed for high performance with minimal maintenance. Made with platinum cured medical grade mineral oil embedded silicone with a Parylene® coating gives the Parasil the unique ability to repel water and lower the co-efficiency of friction. This hydrophobic property is easy to clean and sanitize leaving it odor free
- The Parylene® surface reduces surface friction between the liner/skin and liner/socket which will help in the prevention of skin shear for fragile skin and blister reduction. The Parasil offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Parasil liner's opaque pigment gives a natural "skin tone" to all ethnic skin types

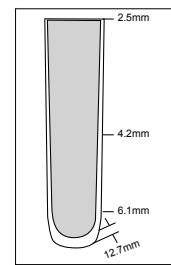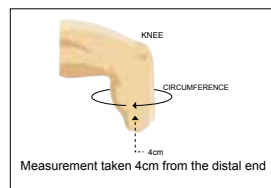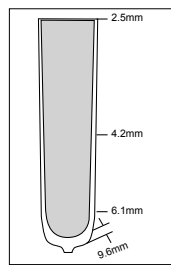

TRANSTIBIAL Locking liner

Product Code	Circumference	Length
3-01-3020L-18	180mm	370mm
3-01-3020L-20	200mm	370mm
3-01-3020L-22	220mm	370mm
3-01-3020L-23.5	235mm	370mm
3-01-3020L-25	250mm	370mm
3-01-3020L-26.5	265mm	370mm
3-01-3020L-28	280mm	370mm
3-01-3020L-30	300mm	370mm
3-01-3020L-32	320mm	370mm
3-01-3020L-34	340mm	370mm
3-01-3020L-36	360mm	370mm

TRANSTIBIAL Cushion liner

Product Code	Circumference	Length
3-01-3020C-18	180mm	360mm
3-01-3020C-20	200mm	360mm
3-01-3020C-22	220mm	360mm
3-01-3020C-23.5	235mm	360mm
3-01-3020C-25	250mm	360mm
3-01-3020C-26.5	265mm	360mm
3-01-3020C-28	280mm	360mm
3-01-3020C-30	300mm	360mm
3-01-3020C-32	320mm	360mm
3-01-3020C-34	340mm	360mm
3-01-3020C-36	360mm	360mm




TRANSFEMORAL Silicone Liner – Parasil (AK)

3-01-3070

3-01-3070

TRANSFEMORAL Silicone Liner (AK)

- Parasil Liner offers a superior durable design for your active patient and is available in both PIN LOCKING and CUSHION
- The silky smooth Parasil Liner is designed for high performance with minimal maintenance. Made with platinum cured medical grade mineral oil embedded silicone with a Parylene® coating gives the Parasil the unique ability to repel water and lower the co-efficiency of friction. This hydrophobic property is easy to clean and sanitize leaving it odor free
- The Parylene® surface reduces surface friction between the liner/skin and liner/socket which will help in the prevention of skin shear for fragile skin and blister reduction. The Parasil offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Parasil liner's opaque pigment gives a natural "skin tone" to all ethnic skin types

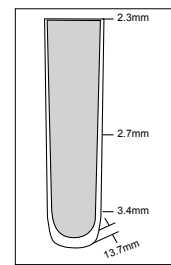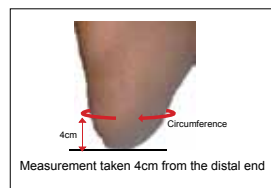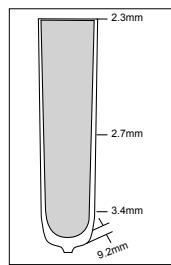

TRANSFEMORAL Locking liner

Product Code	Circumference	Length
3-01-3070L-30	300mm	around 370mm
3-01-3070L-32	320mm	around 370mm
3-01-3070L-34	340mm	around 370mm
3-01-3070L-36	360mm	around 370mm
3-01-3070L-38	380mm	around 370mm
3-01-3070L-40	400mm	around 370mm
3-01-3070L-45	450mm	around 370mm

TRANSFEMORAL Cushion liner

Product Code	Circumference	Length
3-01-3070C-30	300mm	around 370mm
3-01-3070C-32	320mm	around 370mm
3-01-3070C-34	340mm	around 370mm
3-01-3070C-36	360mm	around 370mm
3-01-3070C-38	380mm	around 370mm
3-01-3070C-40	400mm	around 370mm
3-01-3070C-45	450mm	around 370mm

TRANSFEMORAL Silicone Liner – Forte (AK)

3-01-3090

3-01-3090
TRANSFEMORAL Silicone Liner (AK)



- The Forte liner is a silicone gel liner with mineral oil embedded in the silicone providing additional healing factor to the skin. The liner is made of medical grade silicone and durable seamless fabric on the outside for easy donning and doffing application. The liner comes with cushion and locking options. It is suitable for low to moderate activity patients. The matrix which is built into the liner prevents distal elongation of the liner for locking application
- 6 months warranty



3-01-3090L

TRANSFEMORAL Locking liner

Product Code	Circumference	Length
3-01-3090L-30	300mm	around 370mm
3-01-3090L-32	320mm	around 370mm
3-01-3090L-34	340mm	around 370mm
3-01-3090L-36	360mm	around 370mm
3-01-3090L-38	380mm	around 370mm
3-01-3090L-40	400mm	around 370mm
3-01-3090L-45	450mm	around 370mm

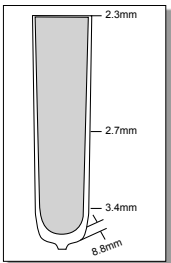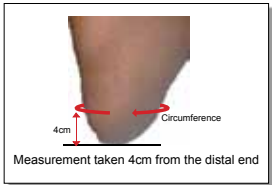


TRANSTIBIAL and TRANSFEMORAL Silicone liner (AK/BK)

3-01-3030

3-01-3030

TRANSTIBIAL and TRANSFEMORAL Silicone liner (AK/BK)

- It made high shore medical grade silicone
- Design for both AK and BK applications
- Fabric on outside for easy don and doff
- 2 types: Locking and Cushioning
- This liner type very suitable for high activities
- Silicone liners has 44cm length totally and it is useful for long stump
- Anti-Bacterial, Anti-Allergic

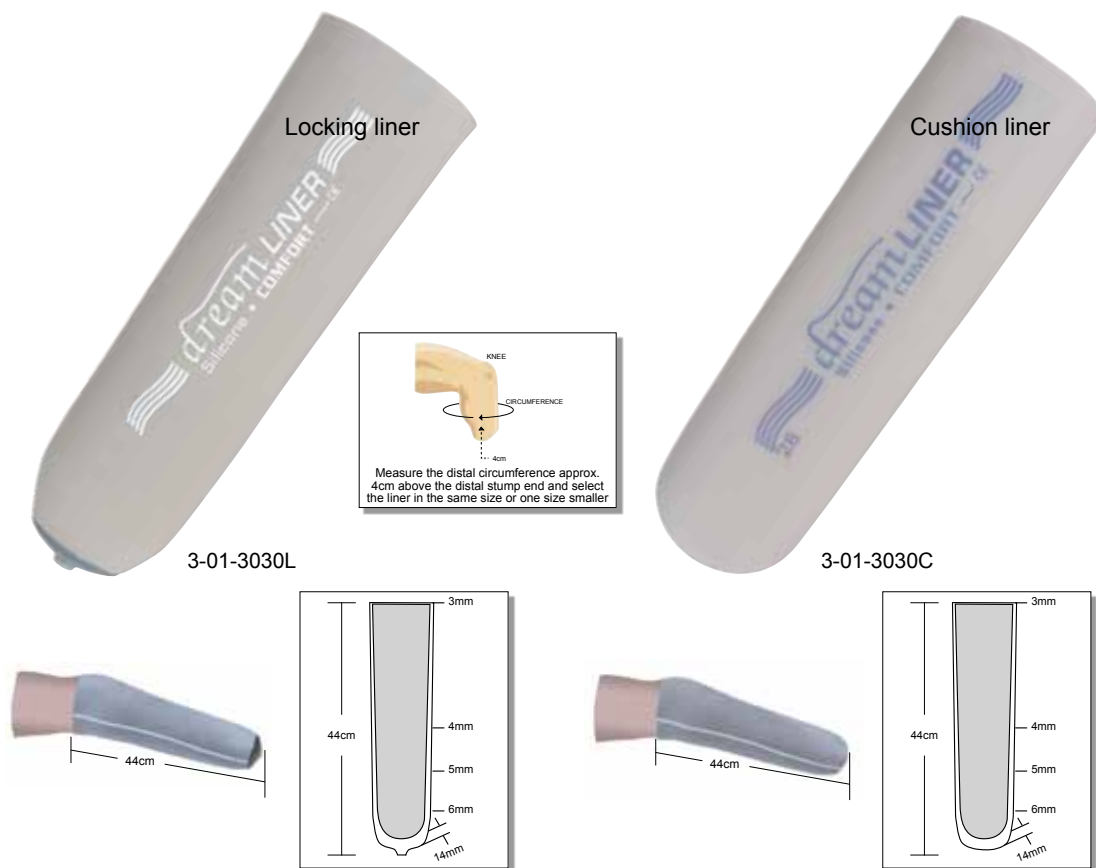

TRANSTIBIAL Locking liner

Product Code	Circum.	Thickness	Length
3-01-3030L-21	210mm	3mm	44cm
3-01-3030L-22	220mm	3mm	44cm
3-01-3030L-23.5	235mm	3mm	44cm
3-01-3030L-25	250mm	3mm	44cm
3-01-3030L-26.5	265mm	3mm	44cm
3-01-3030L-28	280mm	3mm	44cm
3-01-3030L-30	300mm	3mm	44cm
3-01-3030L-32	320mm	3mm	44cm
3-01-3030L-34	340mm	3mm	44cm
3-01-3030L-36	360mm	3mm	44cm
3-01-3030L-38	380mm	3mm	44cm
3-01-3030L-40	400mm	3mm	44cm
3-01-3030L-42	420mm	3mm	44cm
3-01-3030L-44	440mm	3mm	44cm

TRANSTIBIAL Cushion liner

Product Code	Circum.	Thickness	Length
3-01-3030C-21	210mm	3mm	44cm
3-01-3030C-22	220mm	3mm	44cm
3-01-3030C-23.5	235mm	3mm	44cm
3-01-3030C-25	250mm	3mm	44cm
3-01-3030C-26.5	265mm	3mm	44cm
3-01-3030C-28	280mm	3mm	44cm
3-01-3030C-30	300mm	3mm	44cm
3-01-3030C-32	320mm	3mm	44cm
3-01-3030C-34	340mm	3mm	44cm
3-01-3030C-36	360mm	3mm	44cm
3-01-3030C-38	380mm	3mm	44cm
3-01-3030C-40	400mm	3mm	44cm
3-01-3030C-42	420mm	3mm	44cm
3-01-3030C-44	440mm	3mm	44cm

SHIFT clutch lock

3-01-SCL/3-01-SCL-P



3-01-SCL

3-01-SCL-P



New!

3-01-SCL/3-01-SCL-P SHIFT clutch lock

- Smooth gear and bearing mechanism for infinite adjustment
- Can be used for laminated and thermoplastic sockets
- Easily compensate for the variance between pin and socket alignment when use with the 4 hole plate
- Quiet Operation

Product Code	Description	Max. Load
3-01-SCL	Shift Clutch Lock	136kg
3-01-SCL-P	Shift Clutch Lock w/Plate	136kg

Refer to "Instructional video" from our website www.regalprosthesis.com

SHIFT shuttle lock

3-01-SSL/3-01-SSL-P



3-01-SSL

3-01-SSL-P



New!

3-01-SSL/3-01-SSL-P SHIFT shuttle lock

- Serrated pin attachment
- Can be used for laminated and thermoplastic sockets
- Easily compensate for the variance between pin and socket alignment when use with the 4 hole plate

Product Code	Description	Max. Load
3-01-SSL	Shift shuttle Lock	136kg
3-01-SSL-P	Shift shuttle Lock w/Plate	136kg

Refer to "Instructional video" from our website www.regalprosthesis.com

Threaded shuttle lock

3-01-TSL



3-01-TSL

3-01-TSL

Threaded shuttle lock

- Screw in modular design provides easy maintenance
- Cost effective pricing
- Pin length 70mm

Product Code	Description	Max. Load
3-01-TSL	Threaded shuttle lock only	100kg
3-01-TSL-S	Shuttle lock with Stainless 3 prong rotating base	100kg

Lamination clutch lock

3-01-LCL



3-01-LCL

3-01-LCL Lamination clutch lock

- Smooth gear mechanism
- Quite operation
- All-in-one design
- Design for socket lamination
- Cost effective pricing
- Pin length 70mm

Product Code	Description	Max. Load
3-01-LCL	Lamination clutch lock	100kg

Base shuttle lock

3-01-BSL



3-01-BSL-00

3-01-BSL Base shuttle lock

- Made from strong nylon plastic
- Very light design
- Versatile cylindrical body with a 4 hole attachment
- Deep pin well for easy insertion
- Very cost effective



3-01-BSL-01

Product Code	Description	Max. Load
3-01-BSL-00	Base shuttle lock	100kg
3-01-BSL-01	Fabrication Dummy Kit for Shuttle lock	100kg

Refer to "Instructional video" from our website www.regalprosthesis.com

Easy fit

3-01-EF

3-01-EF Easy fit

- Easy to don the AK socket with valve hole
- Ensures good fit between stump and socket
- Easy to use and carry
- Optimum friction reduction between stump and socket
- Durable
- Can be hand-washed



3-01-EF

Product Code	Size	Stump entry circum.
3-01-EF-S	Small	up to 40cm
3-01-EF-M	Medium	41-56cm
3-01-EF-L	Large	57-68cm
3-01-EF-XL	X Large	69-78cm

Screw valve for thermoforming technique

3-01-V8

New!

3-01-V8 Screw valve for thermoforming technique



3-01-V8

Product Code

3-01-V8

XP AK valve

3-01-V7

New!

3-01-V7 XP AK valve

- Simple installation
- Low build height
- For use any activity level
- Durable



3-01-V7

Product Code

3-01-V7

XP valve

3-01-V6



3-01-V6

New!

3-01-V6 XP valve

- Light weight design
- Low build height
- Easy and quick fabrication
- Easy cleaning

Product Code

3-01-V6

Automatic valve

3-01-V5



3-01-V5

New!

3-01-V5 Automatic valve

Product Code

3-01-V5

Prosthetics rubber valve set

3-01-V4



3-01-V4

3-01-V4 Prosthetics rubber valve with seat ring

Product Code

3-01-V4

Flat rubber valve with tube

3-01-V2



3-01-V2

3-01-V2 Flat rubber valve with tube

Product Code

3-01-V2

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Suspension
Systems

Cosmetic, Tools
& Materials

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

Valve kit

3-01-V1



3-01-V1

3-01-V1 Valve kit

Product Code

3-01-V1

AK Suspension sleeve

3-01-AKSS



3-01-AKSS

3-01-AKSS AK-Suspension sleeve

Product Code	Size	Hip size	Socket size	Side
3-01-AKSS-L-M	Medium	66-80cm	40-44cm	Left
3-01-AKSS-R-M	Medium	66-80cm	40-44cm	Right
3-01-AKSS-L-L	Large	76-90cm	44-48cm	Left
3-01-AKSS-R-L	Large	76-90cm	44-48cm	Right
3-01-AKSS-L-XL	X Large	86-100cm	48-52cm	Left
3-01-AKSS-R-XL	X Large	86-100cm	48-52cm	Right
3-01-AKSS-L-XXL	XX Large	96-110cm	52-56cm	Left
3-01-AKSS-R-XXL	XX Large	96-110cm	52-56cm	Right

Suspension sleeve

3-01-KB



3-01-KB

New!

3-01-KB Suspension sleeve

- Suspension sleeve but not for vacuum
- Two way stretch
- Caucasian skin color

Product Code	Size	Length	Below knee	Knee	Above knee
3-01-KB-S	Small	30.5cm	24.5cm	26cm	27.5cm
3-01-KB-M	Medium	32.5cm	28.5cm	30cm	31.5cm
3-01-KB-L	Large	37.5cm	32cm	34.5cm	36.5cm

Eclipse Silicone Suspension Sleeve

3-01-ECSS



3-01-ECSS



3-01-ECSS Eclipse Silicone Suspension Sleeve

The Eclipse Silicone Suspension Sleeve provides extreme durability and comfort for your patient. Imagine not having to replace suspension sleeves every 3 months. A unique combination of materials and smart design has created a superior suspension sleeve that offers high performance and secure fit.

- Reinforced silicone and strategically located Parylene® coating provides maximum durability
- Dual sleeve design with non-slip inner strip of surface provides a maximum airtight seal which eliminates slippage
- This revolutionary patent pending design helps keep the sleeve in place even during the most strenuous activity
- Depending on activity, the outer Eclipse sleeve easily folds down over itself to reduce the length
- The Parylene® coating does not inhibit knee flexion which allows for excellent performance and comfort
- The Eclipse offers a platinum cured medical grade silicone with a Parylene® surface which reduces surface friction between the sleeve and skin helping in the prevention of skin shear for fragile skin and blister reduction
- The smooth and silky coating allows for clothing to naturally hang and not bunch up
- Warranty: 120 day warranty against defect and workmanship

Measure circumference 15cm (6in) above mid-patella and choose the correct size

Product Code	Color	Set contains	Size
3-01-ECSS-S-BK	Black	Long (outer)+Short (inner)	Small (41): 41cm – 49cm
3-01-ECSS-M-BK	Black	Long (outer)+Short (inner)	Medium (49): 49cm – 55cm
3-01-ECSS-L-BK	Black	Long (outer)+Short (inner)	Large (55): 55cm – 62cm
3-01-ECSS-S-TN	Tan	Long (outer)+Short (inner)	Small (41): 41cm – 49cm
3-01-ECSS-M-TN	Tan	Long (outer)+Short (inner)	Medium (49): 49cm – 55cm
3-01-ECSS-L-TN	Tan	Long (outer)+Short (inner)	Large (55): 55cm – 62cm


FAQ & Warranty	Upperlimbs Modular Prosthesis	Orthotics & Rehabilitation	Cosmetic, Tools & Materials	Suspension Systems	Kneejoints & Hip Joints	Feet & Adaptors	Pediatric
----------------	----------------------------------	-------------------------------	--------------------------------	-----------------------	----------------------------	--------------------	-----------

Cosmetic, Tools and Materials



Above knee foam cover

4-01-FA

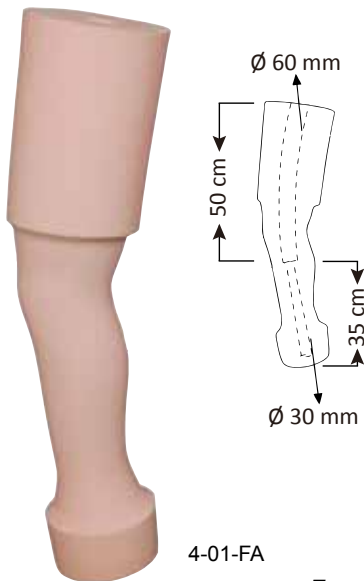

4-01-FA

Above knee foam cover


- Durable foam
- 30° knee flexion
- Preshaped
- Center Hole
- Flame Resistant
- Skin Color

Product Code	Calf Circum.	Side	Hole diameter	Length
4-01-FA-L-36	36cm	Left	60-30mm	85cm
4-01-FA-R-36	36cm	Right	60-30mm	85cm
4-01-FA-L-40	40cm	Left	60-30mm	85cm
4-01-FA-R-40	40cm	Right	60-30mm	85cm
4-01-FA-L-44	44cm	Left	60-30mm	85cm
4-01-FA-R-44	44cm	Right	60-30mm	85cm

Foam covers are available in vacuum pack

Knee disarticulation foam cover

4-01-FKD

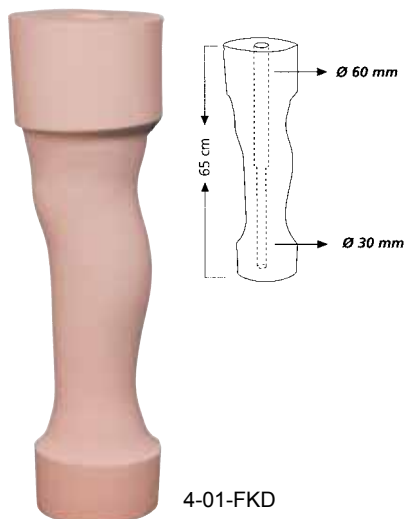

4-01-FKD

Knee disarticulation foam cover

- Durable foam
- Preshaped
- Center Hole
- Flame Resistant
- Skin Color

Product Code	Calf Circum.	Side	Hole diameter	Length
4-01-FKD-L-36	36cm	Left	60-30mm	65cm
4-01-FKD-R-36	36cm	Right	60-30mm	65cm
4-01-FKD-L-40	40cm	Left	60-30mm	65cm
4-01-FKD-R-40	40cm	Right	60-30mm	65cm
4-01-FKD-L-44	44cm	Left	60-30mm	65cm
4-01-FKD-R-44	44cm	Right	60-30mm	65cm

Foam covers are available in vacuum pack

Above knee foam cover (pediatric)

0-01-FAC



0-01-FAC

Above knee foam cover (pediatric)

- Durable foam
- One size
- 30° knee flexion
- Preshaped
- Center Hole
- Flame Resistant
- Skin Color

Product Code	Calf Circum.	Hole diameter	Length
0-01-FAC	30-22cm	22-30mm	60cm


Foam covers are available in vacuum pack

Below knee foam cover

4-01-FB



4-01-FB

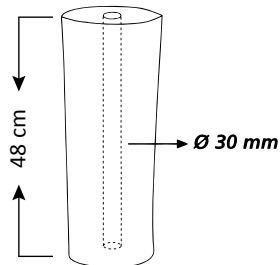

4-01-FB

Below knee foam cover

- Durable foam
- One size
- Unshaped
- Center Hole
- Flame Resistant
- Skin Color



Product Code	Hole diameter	Length
4-01-FB	30mm	48cm


Foam covers are available in vacuum pack

Below knee foam cover (pediatric)

0-01-FBC



0-01-FBC

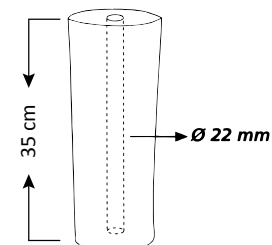

0-01-FBC

Below knee foam cover (pediatric)

- Durable foam
- One size
- Unshaped
- Center Hole
- Flame Resistant
- Skin Color

Product Code	Hole diameter	Length
0-01-FBC	22mm	35cm


Foam covers are available in vacuum pack

Foam cover - Plastazote

4-01-FBP



4-01-FBP

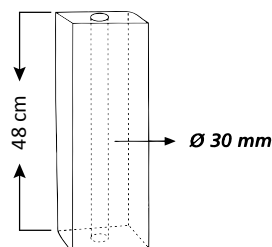

4-01-FBP


Foam cover - Plastazote

- Durable foam
- Unshaped
- Center Hole
- Flame Resistant
- Skin Color
- size:160x160x480 mm

Product Code	Hole diameter	Length
4-01-FBP	30mm	48cm

Cosmetic foam cover for upper limb

4-01-FBE

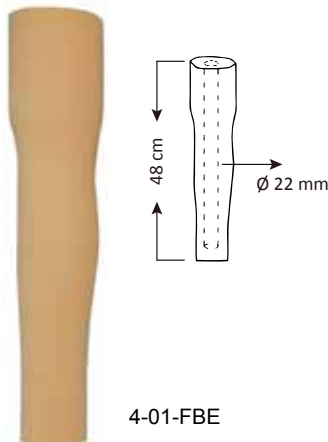

4-01-FBE
Cosmetic foam cover for upper limb
below elbow prosthesis

Product Code	Material	Hole diameter	Length
4-01-FBE	Polyurethane	22mm	47cm

Cosmetic foam covers are available in vacuum pack

Cosmetic foam cover for upper limb

4-01-FSD



4-01-FSD
Cosmetic foam cover for upper limb
shoulder disarticulation upper part

Product Code	Material	Hole diameter	Length
4-01-FSD	Polyurethane	22mm	47+20cm

Cosmetic foam covers are available in vacuum pack

Cosmetic stockings with foot shape

4-01-AKSTO



4-01-AKSTO
Cosmetic stockings with foot shape

- For AK amputation
- One size
- Total length: 70cm

Product Code
4-01-AKSTO

Cosmetic stockings

4-01-BKSTO



4-01-BKSTO
Cosmetic stockings

- For BK amputation
- One size
- Total length: 45cm

Product Code
4-01-BKSTO

Uni-Grip

4-02-UG




Stops Rotation one strip at a time

New!

4-02-UG Uni-Grip

This unique fabric sheet with pressure sensitive adhesive back provides versatile application for prosthetic and orthotic needs. It's biased direction texture provides grip in one direction while providing free glide in the other direction

- Restrict rotation of stump in AK socket
- Suspension for BK and AK socket
- Restrict migration of orthotics
- Open to other ideas...

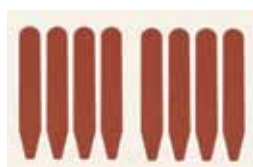

4-02-UG-1



4-02-UG-2

Refer to "Instructional video" from our website www.regalprosthesis.com

Product Code	Description
4-02-UG-1	1 sheet of 8 strips
4-02-UG-2	1 sheet of 2 strips

High quality plastazote

4-02-PZ



4-02-PZ

4-02-PZ

High quality plastazote

- Thermo plastically moldable at 110°C / 230°F on the heating plate, as well as in the convection oven

Product Code	Thickness	Size
4-02-PZ-3	3mm	2m x 1m
4-02-PZ-6	6mm	2m x 1m
4-02-PZ-12	12mm	2m x 1m

High quality polyfoam

4-02-PF



4-02-PF

4-02-PF

High quality polyfoam

- Density 140kg/m3
- Thermo plastically mouldable at 130°C
- Suitable for the production of soft liner sockets

Product Code	Thickness	Size
4-02-PF-2	2mm	2m x 1m
4-02-PF-3	3mm	2m x 1m
4-02-PF-4	4mm	2m x 1m
4-02-PF-5	5mm	2m x 1m
4-02-PF-6	6mm	2m x 1m

Multicolor soft pad


4-02-MP

New!

4-02-MP

Multicolor soft pad

- Thermo plastically moldable at approx. 100-140°C / 212-284°F on the heating plate, and in the oven

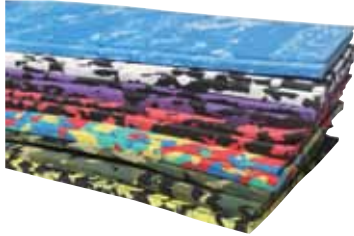


4-02-MP



Product Code	Thickness	Size	Color
4-02-MP-2-WTBK	2mm	2m x 1m	white black
4-02-MP-4-WTBK	4mm	2m x 1m	white black
4-02-MP-2-PUBK	2mm	2m x 1m	purple black
4-02-MP-4-PUBK	4mm	2m x 1m	purple black
4-02-MP-2-MCOL	2mm	2m x 1m	multi color
4-02-MP-4-MCOL	4mm	2m x 1m	multi color
4-02-MP-2-ARMY	2mm	2m x 1m	army
4-02-MP-4-ARMY	4mm	2m x 1m	army
4-02-MP-2-PKBK	2mm	2m x 1m	pink black
4-02-MP-4-PKBK	4mm	2m x 1m	pink black
4-02-MP-2-BLWT	2mm	2m x 1m	blue white
4-02-MP-4-BLWT	4mm	2m x 1m	blue white
4-02-MP-2-YLBK	2mm	2m x 1m	yellow black
4-02-MP-4-YLBK	4mm	2m x 1m	yellow black

Polypropylene (160° -170°)

4-02-PP

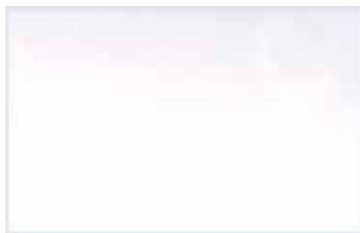

4-02-PP

4-02-PP

Polypropylene (160° -170°)

- For use orthotic applications like AFO, corset, splints

Product Code	Thickness	Size	Color
4-02-PP-2	2mm	2m x 1m	White
4-02-PP-3	3mm	2m x 1m	White
4-02-PP-4	4mm	2m x 1m	White
4-02-PP-5	5mm	2m x 1m	White

Polyethylene (160° -170°)

4-02-PE



4-02-PE

4-02-PE

Polyethylene (160° -170°)

- For use orthotic applications like AFO, corset, splints

Product Code	Thickness	Size	Color
4-02-PE-2	2mm	2m x 1m	White
4-02-PE-3	3mm	2m x 1m	White
4-02-PE-4	4mm	2m x 1m	White
4-02-PE-5	5mm	2m x 1m	White

Check socket material (160°)

4-02-CSM



4-02-CSM

4-02-CSM Check socket material (160°)

- For BK / AK as for test socket

Product Code	Thickness	Size
4-02-CSM-10	10mm	400 x 400mm
4-02-CSM-12	12mm	400 x 400mm
4-02-CSM-15	15mm	400 x 400mm

Elastomer soft socket material (160°)

4-02-ESM



4-02-ESM

4-02-ESM Elastomer soft socket material (160°)

- For AK prosthetics as an inner soft socket

Product Code	Thickness	Size
4-02-ESM-6	6mm	400 x 400mm
4-02-ESM-9	9mm	400 x 400mm
4-02-ESM-12	12mm	400 x 400mm

Orthotic alignment kit

4-02-OAK



4-02-OAK

4-02-OAK Orthotic alignment kit

Product Code

4-02-OAK

Above knee casting brim set (Premium)

4-02-PAKCB



4-02-PAKCB

4-02-PAKCB
Above knee casting brim set
(Premium)

- 22 left & 22 right, total 44 pcs
- For getting measurement from AK residual limb

Product Code

4-02-PAKCB

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Suspension
Systems

Cosmetic, Tools
& Materials

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

Refer to "Instructional video" from our website www.regalprosthesis.com

The Doff N' Donner

4-02-DND



4-02-DND



4-02-DND The Doff N' Donner

- The DoffN' Donner takes the work out of applying compression stockings
- It conforms to the garment, which increases the life of the garment
- Invented with a "push" mentality rather than the "pull" strategy
- The DoffN' Donner easily rolls over bandages, wet limbs and even scar management pads
- Don & Doff Compression garments of all shapes and compression without frustration

Product Code

4-02-DND

POLE for the Doff'N Donner

4-02-PDND



4-02-PDND

New!

4-02-PDND POLE for the Doff'N Donner

- The POLE is the device that "loads" or helps roll the stocking onto the DoffN' Donner
- The base of this pole allows for both foot hold as well as suction
- We would recommend this item for those that do not have a surface to use a suction cup on
- Please note: The power of the suction cup may leave behind a ring when used on a stained finished surface

Product Code

4-02-PDND

CONE for the Doff'N Donner

4-02-CDND



4-02-CDND

New!

4-02-CDND CONE for the Doff'N Donner

- The CONE is the device that "loads" or helps roll the stocking onto the Doff N' Donner
- This device works better for open toe garments, however we have found it to work well with all garment types
- The base of the cone attaches to any clean flat surface using a suction cup
- Please note: The power of the suction cup may leave behind a ring when used on a stained finished surface

Product Code

4-02-CDND

Regal silicone casting resin

4-02-CM



4-02-CM

4-02-CM
Regal silicone casting resin
(for semi-custom made casting only)

- High elasticity
- Good for all types of residual limbs
- 1kg per pail

Product Code

4-02-CM

Regal silicone casting hardener

4-02-CMH



4-02-CMH

4-02-CMH
Regal silicone casting hardener

- High elasticity
- Good for all types of stump
- To be used with Regal casting resin (4-02-CM)
- 100g per bottle

Product Code

4-02-CMH

Regal premium grade color pigments

4-02-PIG



4-02-PIG

4-02-PIG
Regal premium grade color pigments

- Packaging: 1kg / 100g
- For socket lamination
- For color mixing

Product Code	Color
4-02-PIG-CA-1kg	Caucasian
4-02-PIG-CA-100g	Caucasian
4-02-PIG-BW-1kg	Brown
4-02-PIG-BW-100g	Brown
4-02-PIG-RE-1kg	Red
4-02-PIG-RE-100g	Red
4-02-PIG-BL-1kg	Blue
4-02-PIG-BL-100g	Blue
4-02-PIG-YL-1kg	Yellow
4-02-PIG-YL-100g	Yellow
4-02-PIG-WT-1kg	White
4-02-PIG-WT-100g	White

Flexible PU foam (Premium quality)

4-02-PUFA



4-02-PUFA

4-02-PUFA

Flexible PU foam (Premium quality)

- Premium quality
- mixing ratio: 2 PUFA:1 PUFB
- 1kg per pail

Product Code

4-02-PUFA

Flexible PU foam hardener

4-02-PUFB



4-02-PUFB

4-02-PUFB

Flexible PU foam hardener

- Hardener
- Similar to Pedilen Flexible Foam 300
- 1kg per pail

Product Code

4-02-PUFB

Plaster bandages

4-02-PB



4-02-PB

4-02-PB

Plaster bandages

- Length 460cm per roll
- Standard quality
- 10 rolls per pack

Product Code	Width
4-02-PB-75	75mm
4-02-PB-150	150mm

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

NSP braids

4-02-NSP



4-02-NSP

New!

4-02-NSP NSP braids

DYNAMIC STIFFNESS with 10 X THE IMPACT RESISTANCE

NSP is a unique natural fiber that is lighter, stronger, offers ½ the stiffness and 10X the impact resistance than carbon. The stiffness can be regulated by fiber orientation to be equal to carbon. NSP can be stiff or flexible depending upon practitioner preference.

NSP wets out fast because of its high fiber interface and low void content. Translated the resin sticks and air does not get trapped easily.

Unlike carbon, NSP is a green eco-friendly natural fiber that creates no flying fibers and no itch which technicians and practitioners love

- Provides great fatigue strength for the socket
- Tested in socket form under ISO guidelines for patients up to 500 lbs
- Available in 3, 4, 5, 6, 8, 10 & 12 inch width & 10 and 50 feet length
- NSP = WHITE - 4-02-NSP

Product Code	Description	Width	Length
4-02-NSP-3-10	NSP Braid	3 inch	10 feet
4-02-NSP-4-10	NSP Braid	4 inch	10 feet
4-02-NSP-5-10	NSP Braid	5 inch	10 feet
4-02-NSP-6-10	NSP Braid	6 inch	10 feet
4-02-NSP-8-10	NSP Braid	8 inch	10 feet
4-02-NSP-10-10	NSP Braid	10 inch	10 feet
4-02-NSP-12-10	NSP Braid	12 inch	10 feet
4-02-NSP-3-50	NSP Braid	3 inch	50 feet
4-02-NSP-4-50	NSP Braid	4 inch	50 feet
4-02-NSP-5-50	NSP Braid	5 inch	50 feet
4-02-NSP-6-50	NSP Braid	6 inch	50 feet
4-02-NSP-8-50	NSP Braid	8 inch	50 feet
4-02-NSP-10-50	NSP Braid	10 inch	50 feet
4-02-NSP-12-50	NSP Braid	12 inch	50 feet

NSP with
Extra White
Pigment



NSP with
Green
Pigment



NSP with
BLACK
Pigment



NSP with
Tshirt transfer
design



PLEASE NOTE:

To achieve a custom color like the socket color shown above in lime green – color pigment was added to the process. Many different colors can be created by adding this pigment.

The NSP material can create incredible color combinations when using pigment or tshirt transfer design.

Refer to "Instructional video" from our website www.regalprosthesis.com

NSP Uni-directional 3k tape

4-02-NSP-3K



4-02-NSP-3K

4-02-NSP-3K NSP Uni-directional 3k tape

Designed to enhance longitudinal and brim strength and stiffness

NSP Unidirectional tape is a green eco-friendly natural fiber that creates no flying fibers and no itch

- Provides great fatigue strength for the socket
- Green eco-friendly natural fiber
- Is lighter than carbon tape
- Is stronger than carbon tape
- ½ the stiffness of carbon tape
- Offers 10x the impact resistance than carbon tape
- Keeps labs clean due to enhanced fiber saturation
- No itch
- Tested in socket form under ISO guidelines for patients up to 500 lbs

Product Code	Description	Width	Length
4-02-NSP-3K-1-30	NSP Uni-directional 3k tape	1 inch	30 feet
4-02-NSP-3K-3-30	NSP Uni-directional 3k tape	3 inch	10 feet

Carbon braid

4-02-C-TB



4-02-C-TB

4-02-C-TB Carbon braid

- 12K Strand Count
- Designed for improved polyester and fiber interface
- Decreased fraying (less dust)

Product Code	Description	Width	Length
4-02-C-TB-3-10	Carbon Tube Braid	3 inch	10 feet
4-02-C-TB-3-50	Carbon Tube Braid	3 inch	50 feet
4-02-C-TB-4-10	Carbon Tube Braid	4 inch	10 feet
4-02-C-TB-4-50	Carbon Tube Braid	4 inch	50 feet
4-02-C-TB-5-10	Carbon Tube Braid	5 inch	10 feet
4-02-C-TB-5-50	Carbon Tube Braid	5 inch	50 feet
4-02-C-TB-6-10	Carbon Tube Braid	6 inch	10 feet
4-02-C-TB-6-50	Carbon Tube Braid	6 inch	50 feet
4-02-C-TB-8-10	Carbon Tube Braid	8 inch	10 feet
4-02-C-TB-8-50	Carbon Tube Braid	8 inch	50 feet
4-02-C-TB-10-10	Carbon Tube Braid	10 inch	10 feet
4-02-C-TB-10-50	Carbon Tube Braid	10 inch	50 feet

PARALEX carbon + NSP tube braid

4-02-PARALEX-TB



4-02-PARALEX-TB



4-02-PARALEX-TB

PARALEX carbon + NSP tube braid

Combining the world's highest tech fibers...

PARALEX = CARBON + NSP
= Stiffness + Strength
= ½ the weight
= Thin Laminations
= Durability

- The NSP material within the Paralex can create incredible color combinations when using pigment

Cutting PARALEX Fiber:

- Before cutting PARALEX fiber, tape the edge of the fiber with masking tape then cut down the middle of the tape line.
- Use a scissor with a serrated edge.
- If you attempt to cut PARALEX fiber with a smooth bladed scissor, use very sharp scissors then fold the braid over on itself-(doubling up) greatly helps

Lamination:

- For the best edging results, the use of Perlon, nylons or 1/2oz Dacron felt is ideal as the first layer and the last layer.

For best results use 80:20 resin

Grinding:

- Use high speed and low pressure.
- Only use sharp abrasives, dull abrasives will create friction and will "burn off" the resin and leave exposed fibers at your edge

PARALEX Technical Data

Edging Techniques:

(The use of inner and outer finish layers is required. See example layouts below)

Finishing edge instructions:

Finish as you would any composite by wet sanding using 220 grit and finish with 400 grit wet/dry sanding paper.

Example Layout TT:

- Typical TT
- 1 layer of 1/2oz felt or 2 layers of Perlon or nylon
- 1 layer of PARALEX
- 2 layers of Nyglass or equivalent.
- 1 layer of PARALEX
- 2 layers of sheer nylon or equivalent

Example Layout TF:

- Typical TF
- 1 layer of 1/2oz felt or 2 layers of Perlon or nylon
- 1 layer of PARALEX
- 2 layers of Nyglass or equivalent.
- 1 layer of PARALEX
- 2 layers of sheer nylon or equivalent

Product Code	Description	Width	Length
4-02-PARALEX-TB-3-10	PARALEX Carbon + NSP Tube Braid	3 inch	10 feet
4-02-PARALEX-TB-3-50	PARALEX Carbon + NSP Tube Braid	3 inch	50 feet
4-02-PARALEX-TB-4-10	PARALEX Carbon + NSP Tube Braid	4 inch	10 feet
4-02-PARALEX-TB-4-50	PARALEX Carbon + NSP Tube Braid	4 inch	50 feet
4-02-PARALEX-TB-5-10	PARALEX Carbon + NSP Tube Braid	5 inch	10 feet
4-02-PARALEX-TB-5-50	PARALEX Carbon + NSP Tube Braid	5 inch	50 feet
4-02-PARALEX-TB-6-10	PARALEX Carbon + NSP Tube Braid	6 inch	10 feet
4-02-PARALEX-TB-6-50	PARALEX Carbon + NSP Tube Braid	6 inch	50 feet
4-02-PARALEX-TB-8-10	PARALEX Carbon + NSP Tube Braid	8 inch	10 feet
4-02-PARALEX-TB-8-50	PARALEX Carbon + NSP Tube Braid	8 inch	50 feet
4-02-PARALEX-TB-10-10	PARALEX Carbon + NSP Tube Braid	10 inch	10 feet
4-02-PARALEX-TB-10-50	PARALEX Carbon + NSP Tube Braid	10 inch	50 feet

Ultra6 carbon tube braid



4-02-ULTC-TB



4-02-ULTC-TB

Ultra6 carbon tape braid

- Light, Durable & Strong
- Very conforming
- Wets out nicely for a smoother finish
- Comparable cost to standard carbon
- Thinner and lighter than standard carbon without sacrificing strength

Product Code	Description	Width	Length
4-02-ULTC-TB-3-10	Ultra6 carbon tube braid	3 inch	10 feet
4-02-ULTC-TB-3-50	Ultra6 carbon tube braid	3 inch	50 feet
4-02-ULTC-TB-4-10	Ultra6 carbon tube braid	4 inch	10 feet
4-02-ULTC-TB-4-50	Ultra6 carbon tube braid	4 inch	50 feet
4-02-ULTC-TB-5-10	Ultra6 carbon tube braid	5 inch	10 feet
4-02-ULTC-TB-5-50	Ultra6 carbon tube braid	5 inch	50 feet
4-02-ULTC-TB-6-10	Ultra6 carbon tube braid	6 inch	10 feet
4-02-ULTC-TB-6-50	Ultra6 carbon tube braid	6 inch	50 feet
4-02-ULTC-TB-8-10	Ultra6 carbon tube braid	8 inch	10 feet
4-02-ULTC-TB-8-50	Ultra6 carbon tube braid	8 inch	50 feet
4-02-ULTC-TB-10-10	Ultra6 carbon tube braid	10 inch	10 feet
4-02-ULTC-TB-10-50	Ultra6 carbon tube braid	10 inch	50 feet

Ultra6 carbon tape

4-02-ULTC-TP



4-02-ULTC-TP

4-02-ULTC-TP

Ultra6 carbon tape

- New lighter weight unidirectional carbon tape for optimum strength

Product Code	Description	Width	Length
4-02-ULTC-TP-1-30	Ultra6 Carbon Tape	1 inch	30 feet
4-02-ULTC-TP-3-30	Ultra6 Carbon Tape	3 inch	30 feet
4-02-ULTC-TP-3-100	Ultra6 Carbon Tape	3 inch	100 feet

Cotton stockinette

4-02-STO-COT



4-02-STO-COT

4-02-STO-COT Cotton stockinette

- For reinforcement of laminating resins

Product Code	Width	Weight	Length
4-02-STO-COT-060	60mm	1kg / roll	~29m / roll
4-02-STO-COT-090	90mm	1kg / roll	~20m / roll
4-02-STO-COT-130	130mm	1kg / roll	~14m / roll
4-02-STO-COT-150	150mm	1kg / roll	~12m / roll
4-02-STO-COT-250	250mm	1kg / roll	~7m / roll
4-02-STO-COT-300	300mm	1kg / roll	~6m / roll

Perlon stockinette

4-02-STO-PER



4-02-STO-PER

4-02-STO-PER Perlon stockinette

- For reinforcement of laminating resins
- Very elastic

Product Code	Width	Weight	Length
4-02-STO-PER-080	80mm	1kg / roll	~52m / roll
4-02-STO-PER-100	100mm	1kg / roll	~44m / roll
4-02-STO-PER-120	120mm	1kg / roll	~35m / roll
4-02-STO-PER-130	130mm	1kg / roll	~34m / roll
4-02-STO-PER-150	150mm	1kg / roll	~30m / roll
4-02-STO-PER-180	180mm	1kg / roll	~26m / roll
4-02-STO-PER-250	250mm	1kg / roll	~20m / roll

Fiberglass stockinette

4-02-STO-FG



4-02-STO-FG

4-02-STO-FG Fiberglass stockinette

- For reinforcement of laminating resins
- Highly elastic
- Blend of fiberglass

Product Code	Width	Weight	Length
4-02-STO-FG-080	80mm	1kg / roll	~48m / roll
4-02-STO-FG-100	100mm	1kg / roll	~35m / roll
4-02-STO-FG-130	130mm	1kg / roll	~30m / roll
4-02-STO-FG-180	180mm	1kg / roll	~20m / roll
4-02-STO-FG-250	250mm	1kg / roll	~15m / roll

Premium quality PVA film

4-02-PVA-P



4-02-PVA-P

4-02-PVA-P Premium quality PVA film

Product Code	Width	Thickness	Length
4-02-PVA-P100	1000mm	0.05mm	100m

Economy PVA film

4-02-PVA-E




4-02-PVA-E

4-02-PVA-E Economy PVA film

Product Code	Width	Thickness	Length
4-02-PVA-E50	1114mm	0.08mm	50m
4-02-PVA-E100	1114mm	0.08mm	100m

PVA bag

4-02-PVAB

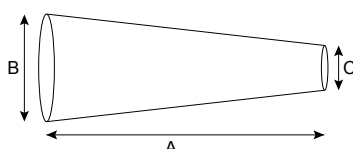

4-02-PVAB

4-02-PVAB PVA bag

- With conical shape for acrylic and polyester laminating resins
- Thickness : 0.08mm
- Water soluble
- 10pcs per pack

Product Code	Dimension A	Dimension B	Dimension C
4-02-PVAB-XS	600mm	130mm	40mm
4-02-PVAB-S	700mm	190mm	50mm
4-02-PVAB-M	700mm	270mm	50mm
4-02-PVAB-L	1000mm	190mm	50mm
4-02-PVAB-XL	1000mm	260mm	50mm

Regal Silicone glue

4-02-SG

Pediatric



4-02-SG-85



4-02-SG-14.2

4-02-SG
Regal Silicone glue

- one component adhesive made specifically for bonding RTV silicone rubber to silicone rubber and other materials
- can be used for acrylic nails

Product Code	Weight
4-02-SG-14.2	14.2gm
4-02-SG-85	85gm

Regal Silicone putty

4-02-SP

Kneejoints & Hip Joints



4-02-SP

4-02-SP
Regal Silicone putty

- It uses for casting purpose



4-02-SP-A


4-02-SP-B

Product Code	Description	Net Contents
4-02-SP-A	Silicone putty (Base)	500 ml
4-02-SP-B	Silicone putty (Catalyst)	500 ml

Deep drawing tool set


4-02-DDS

Orthotics & Rehabilitation

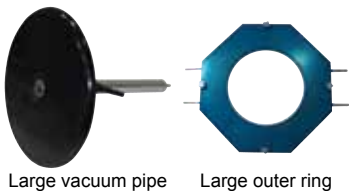

Large vacuum pipe

Large outer ring

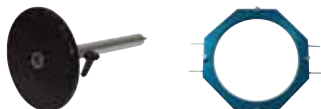

Small vacuum pipe

Small outer ring

4-02-DDS



Vacuum forming tool kit stand

4-02-DDS
Deep drawing tool set

Product Code
4-02-DDS

Upperlimbs
Modular Prosthesis

FAQ & Warranty

Connection cap for Regal single axis foot

1-01-FT-01



1-01-FT-01

1-01-FT-01 Connection cap for Regal single axis foot

Product Code	Size
1-01-FT-01-21	21cm
1-01-FT-01-22	22cm
1-01-FT-01-23	23cm
1-01-FT-01-24	24cm
1-01-FT-01-25	25cm
1-01-FT-01-26	26cm
1-01-FT-01-27	27cm
1-01-FT-01-28	28cm

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Connection plate for Regal SACH foot

1-01-FTT-01



1-01-FTT-01

1-01-FTT-01 Connection plate for Regal SACH foot

- Glue to Regal cosmetic foam cover


Product Code
1-01-FTT-01

Suspension
Systems

Cosmetic, Tools
& Materials

Grease pencil

4-02-GP

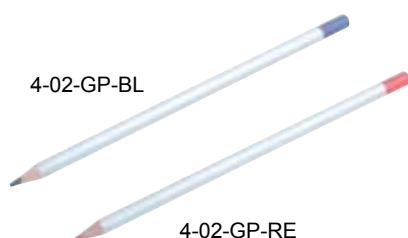

4-02-GP-BL

4-02-GP-RE

4-02-GP Grease pencil

- For marks on wet plaster molds and plaster casts
- 10 pcs per pack

Product Code	Color
4-02-GP-BL	Blue
4-02-GP-RE	Red

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials


Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Plastic loop

4-02-PL

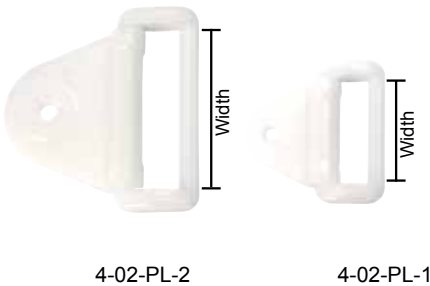


4-02-PL
Plastic loop

- 50 pcs per box

Product Code	Size	Width (inner loop)
4-02-PL-1	Small	2.5cm
4-02-PL-2	Large	3.8cm

Plastic bolt

4-02-PBOLT

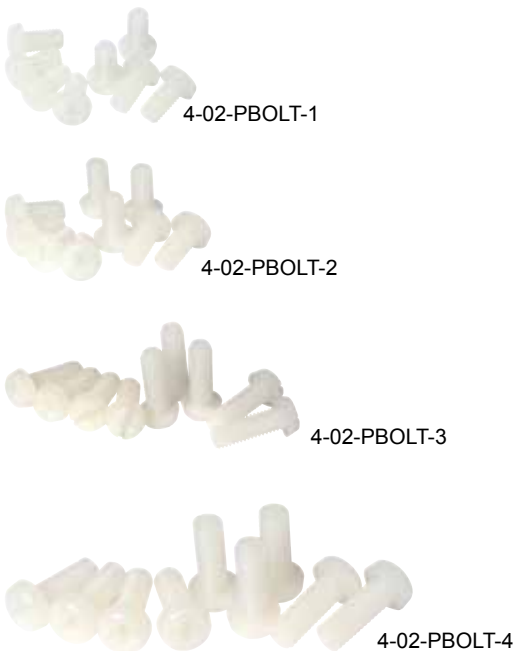

4-02-PBOLT
Plastic bolt

- 50 pcs per box

Product Code	Size
4-02-PBOLT-1	M5 x 15
4-02-PBOLT-2	M5 x 20
4-02-PBOLT-3	M6 x 16
4-02-PBOLT-4	M6 x 20

Pipe cutter

4-02-PC



4-02-PC
Pipe cutter

- cut tube size (3-32mm)
- Suitable for stainless steel, aluminium and copper tube

4-02-PC

Product Code
4-02-PC

Wooden spatula

4-02-WS



4-02-WS-150

4-02-WS-160

4-02-WS Wooden spatula

- Beech wood
- 50 pcs per pack

Product Code	Length	Width	Thickness
4-02-WS-150	150mm	18mm	1.6mm
4-02-WS-160	160mm	10mm	2mm

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Thread locker

4-02-LOC242/4-02-LOC271



4-02-LOC242

4-02-LOC271

4-02-LOC242/4-02-LOC271 Thread locker

- Designed for precision metal fasteners
- Protect threads from rust and corrosion

Product Code	Net Contents	Strength
4-02-LOC242	50ml / bott	Medium
4-02-LOC271	50ml / bott	Strong

Suspension
Systems

Cosmetic, Tools
& Materials

Medical adhesive

4-02-MAD



4-02-MAD

4-02-MAD Medical adhesive

- Very good medical grade adhesive for holding silicone prostheses onto the stump

Product Code

4-02-MAD

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

Sanding drum

4-02-SD



4-02-SD01



4-02-SD02M / 4-02-SD02I



4-02-SD03M / 4-02-SD03I

4-02-SD Sanding drum

- Sanding tools for routers

Product Code	Size	Thread
4-02-SD01	Φ72.5 x L200mm	-
4-02-SD02M	Φ45 x L70mm	M16
4-02-SD02I	Φ45 x L70mm	5/8 inch
4-02-SD03M	Φ25 x L45mm	M16
4-02-SD03I	Φ25 x L45mm	5/8 inch

Sanding sleeve

4-02-SL



4-02-SL01



4-02-SL02



4-02-SL03

4-02-SL Sanding sleeve

- Sanding tools for routers
- For 4-02-SD use
- When ordering, please indicate grit 40#, 60#, 80# or 100#

Product Code	Size	Grit #	Compatible with
4-02-SL01-040	Φ72.5 x L200mm	40	4-02-SD01
4-02-SL01-060	Φ72.5 x L200mm	60	4-02-SD01
4-02-SL01-080	Φ72.5 x L200mm	80	4-02-SD01
4-02-SL01-100	Φ72.5 x L200mm	100	4-02-SD01
4-02-SL02-040	Φ45 x L70mm	40	4-02-SD02
4-02-SL02-060	Φ45 x L70mm	60	4-02-SD02
4-02-SL02-080	Φ45 x L70mm	80	4-02-SD02
4-02-SL02-100	Φ45 x L70mm	100	4-02-SD02
4-02-SL03-060	Φ25 x L45mm	60	4-02-SD03
4-02-SL03-080	Φ25 x L45mm	80	4-02-SD03
4-02-SL03-100	Φ25 x L45mm	100	4-02-SD03

Friction wrist connector (for prosthetic socket lamination)

4-02-WC

4-02-WC

Friction wrist connector for prosthetic socket lamination

- Constant friction applied with a unique friction adjustment collar
- Consists of:
 - a friction unit connector (choice of M12x1.5 or 1/2"-20)
 - a friction adjustment collar
 - a thread stud (choice of M12x1.5 or 1/2"-20)



Product Code	Shape	Thread	Stud	Size
4-02-WC-MT-L	Oval shape	M12x1.5	Male/Female	L (550x395x160mm)
4-02-WC-IM-L	Oval shape	1/2"-20	Male/Female	L (550x395x160mm)
4-02-WC-MT-M	Oval shape	M12x1.5	Male/Female	M (490x380x160mm)
4-02-WC-IM-M	Oval shape	1/2"-20	Male/Female	M (490x380x160mm)
4-02-WC-MT-S	Oval shape	M12x1.5	Male/Female	S (380x280x160mm)
4-02-WC-IM-S	Oval shape	1/2"-20	Male/Female	S (380x280x160mm)

Hand Wrist unit with constant friction

4-02-WH

4-02-WH

Hand wrist unit with constant friction

- Internal thread: M12 x 1.5



Product Code	Size (Outer diameter (Ø))
4-02-WH-40	40mm
4-02-WH-45	45mm
4-02-WH-50	50mm

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Plaster blade handle

4-02-PBHE



4-02-PBHE

4-02-PBHE
Plaster blade handle

Product Code
4-02-PBHE

Plaster blade

4-02-PBFL



4-02-PBFL

4-02-PBFL
Plaster blade

Product Code	Blade Shape	Width
4-02-PBFL	Flat	40mm

Plaster blade

4-02-PBHR



4-02-PBHR

4-02-PBHR
Plaster blade

Product Code	Blade Shape	Width
4-02-PBHR	Half round	40mm

Plaster blade

4-02-PBRD



4-02-PBRD

4-02-PBRD
Plaster blade

Product Code	Blade Shape	Diameter
4-02-PBRD	Round	150mm

Orthotics and Rehabilitation



Posterior leaf spring AFO

5-01-AFO001



5-01-AFO001

5-01-AFO001 Posterior leaf spring AFO

- The Posterior Leaf Spring AFO fulfills immediate patient care needs, assists with drop foot control
- Provides medial/lateral stability at ankle and can be utilized as an evaluation tool
- This semi-rigid device is manufactured with orthopedic grade polypropylene and a 1½ hook and loop proximal calf strap with a closed cell foam shin pad. Right or left needs to be specified

Product Code	Size	Side	Foot Length
5-01-AFO001-L-S	Small	Left	23cm
5-01-AFO001-R-S	Small	Right	23cm
5-01-AFO001-L-M	Medium	Left	24.5cm
5-01-AFO001-R-M	Medium	Right	24.5cm
5-01-AFO001-L-L	Large	Left	26.5cm
5-01-AFO001-R-L	Large	Right	26.5cm
5-01-AFO001-L-XL	X Large	Left	28cm
5-01-AFO001-R-XL	X Large	Right	28cm

Swedish AFO

5-01-AFO002



5-01-AFO002

5-01-AFO002 Swedish AFO

- The Swedish AFO fulfills immediate patient care needs, assists with drop foot control
- Provides moderate spring assist and can be utilized as an evaluation tool
- This semi-rigid device is manufactured with orthopedic grade polypropylene and a 1½ hook and loop proximal calf strap with a closed cell foam shin pad. Right or left needs to be specified

Product Code	Size	Side	Foot Length
5-01-AFO002-L-S	Small	Left	23cm
5-01-AFO002-R-S	Small	Right	23cm
5-01-AFO002-L-M	Medium	Left	24.5cm
5-01-AFO002-R-M	Medium	Right	24.5cm
5-01-AFO002-L-L	Large	Left	26.5cm
5-01-AFO002-R-L	Large	Right	26.5cm
5-01-AFO002-L-XL	X Large	Left	28cm
5-01-AFO002-R-XL	X Large	Right	28cm

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Ridge Flex Ankle Joint

5-01-RF

5-01-RF Ridge Flex Ankle Joint



5-01-RFD100

- The RidgeFlex Ankle Joint represents optimal balance between flexibility and rigidity due to the utilization of a proprietary material called Prepalon
- Prepalon is a continuous fiber and elastomer blend technology (CFET)
- The use of flexible elastomer provides optimum tear strength, preventing fatigue failure
- Continuous fiber blended within the elastomer provides structure and rigidity, preventing the joint from stretching
- This technology eliminates the need for strings which do not provide any support in rotation or A-P movement
- Material technology allows for bold concepts which led to the creation of the tendon design. The tendon design offers a smooth movement while providing proper support in the sagittal plane and rotational control
- Combination of flexible elastomer and continuous fiber blend for rigidity allows for optimum balance of smooth movement and rigid support
- No need to compromise the tear strength by using a higher density elastomer for more support. Higher density in elastomer means weaker tear strength, thus compromising the structural integrity of the joint



5-01-RFS100



CFET

continuous fiber and elastomer blend technology

Product Code	Description
5-01-RFD100	Dorsiflexion Assist
5-01-RFS100	Standard
5-01-RFS100DM	Dummy – (Dorsiflexion Assist & Standard)

Knight armor - Nightsplint

5-01-PFSS

5-01-PFSS

Knight armor - Nightsplint

- The Knight Armor provides relief from pain and discomfort of plantar fasciitis and achilles tendonitis
- Applies consistent stretch to plantar fascia by holding foot in desired dorsi-flexion angle
- Dorsi-flexion assist/resist straps, easy push lock buckles, padded forefoot strap, padded proximal calf strap, non-slip walking pad, washable soft-foam liner, light weight design for patient comfort, universal for right or left applications



5-01-PFSS

Product Code	Size
5-01-PFSS-S	Small
5-01-PFSS-M	Medium
5-01-PFSS-L	Large

Universal hip joint

5-01-HJ

5-01-HJ

Universal hip joint

- The universal hip joint has a sleek low profile design and can be used for left or right applications
- The uprights are manufactured utilizing electroplated aircraft aluminum and the midsection is machined out of high quality stainless steel
- The joint design allows for simple adjustment in adduction/abduction, flexion/extension, range of motion and length/height
- The adduction/abduction is simply adjusted with hex screw in 5° intervals, set and locked with stainless steel teeth assuring strength and durability
- Flexion/extension is adjusted in 15° intervals. Range of motion is also set in 15° intervals and can be completely locked out at any position or set at any desired range motion within the 15° intervals
- The joint also features slotted uprights to easily adjust length/height. The universal hip joint is also cost effective and can be utilized in a wide variety of orthotic applications



5-01-HJ

Product Code

5-01-HJ

Ratcheting buckle

5-01-RTB

5-01-RTB

Ratcheting buckle

- Ratcheting buckle can be used as straps in any orthotic application
- The ratchet design makes it easy for the patient to don the brace with little effort and the quick release enables the patient to doff the brace easily as well
- Two colors are available.



5-01-RTB

Product Code	Color
5-01-RTB-WT	White
5-01-RTB-BK	Black

Criss-cross hyperextension spinal orthosis

5-01-RB01



5-01-RB01

5-01-RB01

Criss-cross hyperextension spinal orthosis

- The Criss-Cross hyperextension spinal orthosis reduces risk of additional injury from osteoporosis and compression fractures
- This orthosis is for patients needing moderate control
- The Criss-Cross frame is manufactured from light weight polished aluminum and features an additional adjustment hole for increased height/length and medial/lateral positioning
- This orthosis also features a padded posterior pad with a hook and loop closure system and a slim line thoracic pad for additional comfort and patient compliance
- The Criss-Cross is also a high quality and cost effective product
- Stable, spontaneous or traumatic compression fractures, reducing further injury from osteoporosis or compression during therapy and rehab process

Product Code

5-01-RB01

Spring hinged hyperextension TLSO

5-01-SB01



5-01-SB01

5-01-SB01

Spring hinged hyperextension TLSO

- The spring hinged hyperextension TLSO features a 3-point pressure control to dynamically maintain leverage for therapeutic immobilization reducing the load on the anterior vertebral bodies
- The single lumbar strap with ratchet buckle creates a snag free closure assisting in patient compliance and ease of donning and doffing. The overall design reduces kyphotic posture associated with compression fractures
- The dynamic extension is accomplished through the spring-hinged design, coupled, with the ratchet locking mechanism creating the 3-point pressure control
- The orthosis also features a lumbar pad, slim line thoracic pad, adjustment for height and M/L and corrosion resistant aluminum frame nicely polished and finished

Product Code	Size
5-01-SB01-S	Small
5-01-SB01-M	Medium
5-01-SB01-L	Large
5-01-SB01-XL	X Large

Meridian ankle joint & fabrication kit

5-01-MAJ001/5-01-MAJKIT



5-01-MAJ001
Meridian Ankle Joint

5-01-MAJ001/5-01-MAJKIT Meridian ankle joint & fabrication kit

- The Meridian ankle joint is a durable stainless steel ankle joint
- Its sleek and low profile design allows for a light weight and durable ankle joint, which can be set in infinite degrees of motion, without the inconvenience of changing parts
- This joint allows precise degree fine tuning and has the ability to lock in any degree of plantar/dorsi-flexion
- Some uses for the Meridian ankle joint include articulating AFO's, contracture management, solid ankle AFO's, posterior entry AFO's, tone reducing AFO's and floor reaction AFO's



5-01-MAJKIT
Fabrication Kit

Product Code	Description	Packing
5-01-MAJ001-P	Pediatric	Pair
5-01-MAJ001-A	Adult	Pair
5-01-MAJ001-HD	Heavy Duty	Pair
5-01-MAJKIT-P	Pediatric	Kit
5-01-MAJKIT-A	Adult	Kit
5-01-MAJKIT-HD	Heavy Duty	Kit

Variable, Varus, Valgus (V-3) correction joint

5-01-V3CJ



5-01-V3CJ
Variable joint, Adult

5-01-V3CJ Variable, Varus, Valgus (V-3) correction joint

- The V-3 joint has some significant benefits and features which include
- Ease of fabrication with the utilization of the polyurethane fabrication molding dummies
- Allows for complete vertical alignment in the coronal plane, simple varus/valgus correction adjustability, a self-aligning midsection that provides a smooth free motion articulation at any angle, superior rotation stability of AFO
- All of these features provide for an easily adjustable ankle joint system during the fitting process
- In addition, interlocking proximal grooves ensure secure joint placement within plastic
- The V-3 is a durable and light weight design machined from high quality stainless steel
- This joint is also extremely cost effective and saves fabrication and fitting time



5-01-V3CJ-01
Variable Joint Fabrication Molding Dummy

Product Code	Description	Packing
5-01-V3CJ	Adult	Pair
5-01-V3CJ-01	Dummy	Piece

PAC Cervical collar

5-01-COL001



5-01-COL001

5-01-COL001 PAC Cervical collar

- Made of light weight closed cell foam material
- Molded and shaped foam helps conform to any patient
- Provides rigid immobilization to cervical spine
- Hook and loop closure for easy adjustment
- Perforated hole for ventilation

Product Code	Size	Circumference	Height
5-01-COL001-2-S	Small	10-13 inch	2.25 inch
5-01-COL001-2-M	Medium	13-16 inch	2.25 inch
5-01-COL001-2-L	Large	16-19 inch	2.25 inch
5-01-COL001-2-XL	X-Large	19-22 inch	2.25 inch
5-01-COL001-3-S	Small	10-13 inch	3.25 inch
5-01-COL001-3-M	Medium	13-16 inch	3.25 inch
5-01-COL001-3-L	Large	16-19 inch	3.25 inch
5-01-COL001-3-XL	X-Large	19-22 inch	3.25 inch
5-01-COL001-4-S	Small	10-13 inch	4.25 inch
5-01-COL001-4-M	Medium	13-16 inch	4.25 inch
5-01-COL001-4-L	Large	16-19 inch	4.25 inch
5-01-COL001-4-XL	X-Large	19-22 inch	4.25 inch
5-01-COL001-5-S	Small	10-13 inch	5.25 inch
5-01-COL001-5-M	Medium	13-16 inch	5.25 inch
5-01-COL001-5-L	Large	16-19 inch	5.25 inch
5-01-COL001-5-XL	X-Large	19-22 inch	5.25 inch

Pediatric

Feet & Adaptors

Kneejoints & Hip Joints

Suspension Systems

Cosmetic, Tools & Materials

Orthotics & Rehabilitation

Upperlimbs Modular Prosthesis

FAQ & Warranty

Ring lock system - Straight

5-01-RLOC



5-01-RLOC

Product Code	Length	Width	Thickness	Weight	Description
5-01-RLOC-SS	900mm	19mm	4mm	1300g	Stainless steel upright bar
5-01-RLOC-AL	900mm	19mm	6mm	850g	Aluminum upright bar

Swiss lock system - Straight

5-01-SLOC



5-01-SLOC

Product Code	Length	Width	Thickness	Weight	Description
5-01-SLOC-SS	850mm	19mm	4mm	1400g	Stainless steel upright bar
5-01-SLOC-AL	850mm	19mm	6mm	900g	Aluminum upright bar

Offset orthotic joint

5-01-OFFS-SS



5-01-OFFS-SS

Product Code	Length	Width	Thickness	Weight	Description
5-01-OFFS-SS	850mm	19mm	4mm	1300g	Offset orthotic joint (stainless steel)

Ring lock system - Straight (pediatric)

0-01-RLOCC-AL



0-01-RLOCC-AL

Product Code	Length	Width	Thickness	Weight	Description
0-01-RLOCC-AL	740mm	12mm	5mm	400g	Aluminum upright bar

Ring lock system - Straight (pediatric)

0-01-RLOCC-SS



0-01-RLOCC-SS

Product Code	Length	Width	Thickness	Weight	Description
0-01-RLOCC-SS	460mm	12mm	3mm	350g	Stainless steel upright bar

Offset orthotic joint (pediatric)

0-01-OFFSC-SS



0-01-OFFSC-SS

Product Code	Length	Width	Thickness	Weight	Description
0-01-OFFSC-SS	850mm	16mm	4mm	1100g	Offset orthotic joint (stainless steel)

Full length silicone insoles

5-01-FL

5-01-FL

Full length silicone insoles

- Made of medical level silicon rubber
- The full length provides all-over pressure relief over a large surface
- Ideal for plantar pain, diabetic feet, heel pain/spurs, metatarsal, arthritis, post-surgical feet and also in sports activities



5-01-FL

Product Code	Size
5-01-FL-XS	XS
5-01-FL-S	S
5-01-FL-M	M
5-01-FL-L	L
5-01-FL-XL	XL

3/4 Thin pad silicone insoles

5-01-QL

5-01-QL

3/4 Thin pad silicone insoles

- Made of medical level silicon rubber
- Ideal for dress shoes or pumps
- The metatarsal support helps the realignment of metatarsal heads, relieves serve forefoot discomfort or callus buildup
- The position of the softer and blue chips helps relieve the pain of metatarsal



5-01-QL

Product Code	Size
5-01-QL-S	S
5-01-QL-M	M
5-01-QL-L	L

Silicone heel cups for spurs - Central

5-01-HCC

5-01-HCC

Silicone heel cups for spurs - Central

- Made of medical level silicone rubber
- Superior shock absorption to relieve and prevent heel, joint and back pain
- The soft blue dot in the center of the heel provides additional cushioning
- Comfortable and long lasting



5-01-HCC

Product Code	Size
5-01-HCC-S	S
5-01-HCC-M	M
5-01-HCC-L	L

Silicone toe spreader

5-01-TS

5-01-TS

Silicone toe spreader



5-01-TS

Product Code	Size
5-01-TS-S	S
5-01-TS-M	M
5-01-TS-L	L

Upperlimbs Modular Prosthesis



Shell locking elbow joint passive device

6-01-EJ1




6-01-EJ1-CA

6-01-EJ1-DKBW

6-01-EJ1

Shell locking elbow joint passive device


- With passive elbow lock, and upper arm rotation joint

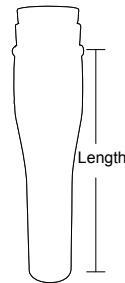

Product Code	Length	Proximal Circumference	Color	Wrist outer circumference	Wrist inner circumference
6-01-EJ1-CA	approx. 287mm	200mm	Caucasian	approx. 54mm	approx. 49mm
6-01-EJ1-DKBW	approx. 287mm	200mm	Dark brown	approx. 54mm	approx. 49mm

Cosmetic prostheses for forearm

6-01-EJ2

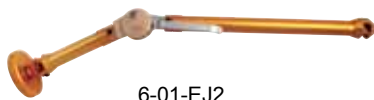

6-01-EJ2

6-01-EJ2

Cosmetic prostheses for forearm

- For above-elbow residual limb, with passive elbow lock, upper arm and forearm rotation, separate friction

Product Code	Side	Threaded Stud
6-01-EJ2-L	Left	M12 x 1.5
6-01-EJ2-R	Right	M12 x 1.5

Cosmetic prostheses for shoulder joint

6-01-SJ1



6-01-SJ1

6-01-SJ1

Cosmetic prostheses for shoulder joint

- For shoulder disarticulation, with passive elbow lock

Product Code	Side	Threaded Stud
6-01-SJ1-L	Left	M12 x 1.5
6-01-SJ1-R	Right	M12 x 1.5


Mechanical hand - voluntary opening (with PVCIH)

6-02-MH1

6-02-MH1

Mechanical hand voluntary opening (with PVCIH)

- With perlon cable through the back of the hand (outer pull)
- With threaded stud and inner hand

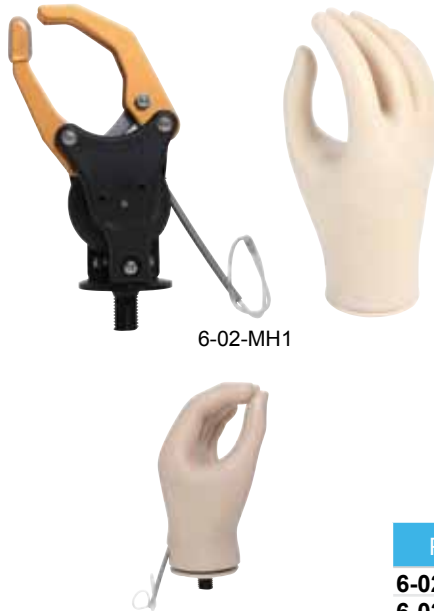

6-02-MH1

Product Code	Side	Size	Threaded stud
6-02-MH1-MT-L-S	Left	S (7")	M12 x 1.5
6-02-MH1-MT-R-S	Right	S (7")	M12 x 1.5
6-02-MH1-MT-L-M	Left	M (7 ¼")	M12 x 1.5
6-02-MH1-MT-R-M	Right	M (7 ¼")	M12 x 1.5
6-02-MH1-MT-L-L	Left	L (7 ¾")	M12 x 1.5
6-02-MH1-MT-R-L	Right	L (7 ¾")	M12 x 1.5
6-02-MH1-MT-L-XL	Left	XL (8" - 8 ¼")	M12 x 1.5
6-02-MH1-MT-R-XL	Right	XL (8" - 8 ¼")	M12 x 1.5
6-02-MH1-IM-L-S	Left	S (7")	1/2"-20
6-02-MH1-IM-R-S	Right	S (7")	1/2"-20
6-02-MH1-IM-L-M	Left	M (7 ¼")	1/2"-20
6-02-MH1-IM-R-M	Right	M (7 ¼")	1/2"-20
6-02-MH1-IM-L-L	Left	L (7 ¾")	1/2"-20
6-02-MH1-IM-R-L	Right	L (7 ¾")	1/2"-20
6-02-MH1-IM-L-XL	Left	XL (8" - 8 ¼")	1/2"-20
6-02-MH1-IM-R-XL	Right	XL (8" - 8 ¼")	1/2"-20

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Suspension
Systems

Cosmetic, Tools
& Materials

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis

FAQ & Warranty

FAQ & Warranty	
Upperlimbs Modular Prosthesis	
Orthotics & Rehabilitation	
Cosmetic, Tools & Materials	
Suspension Systems	
Kneejoints & Hip Joints	
Feet & Adaptors	
Pediatric	

Pediatric

Feet &
Adaptors

Kneejoints &
Hip Joints

Suspension
Systems

Cosmetic, Tools
& Materials

Orthotics &
Rehabilitation

Upperlimbs
Modular Prosthesis



REGAL Prosthesis Ltd.

Room 3D, Tower F, Mai Luen Industrial Bldg.,
23-31 Kung Yip Street, Kwai Chung, NT, Hong Kong.
T : +852 2771-8991 F : +852 2384-5948



www.regalprosthesis.com



info@regalprosthesis.com



Regal Prosthesis Limited



Regal Prosthesis Ltd

